company profile

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RESISTA, trademark owned by Ing. Carlo Alberto Issoglio & C. S.r.l. , identifies the full range of products designed by the company to meet the needs of dentists and dental technicians.

RESISTA is known all over the world since 1946, and synonymous of guarantee, balance and continuity.

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The company has one of the most up-to-date manufacturing site for the medical devices, located in Omegna, a beautiful town on Orta Lake, Italy.

With commitment, merit and pride, the technical and scientific staff of the R&D department is daily involved to promote innovation and the improvement of our products and processes.

Over the years, the **RESISTA** brand has also become synonymous of fast and qualitative response to market changes and the needs of modern dentistry.

Quality, Price and Service summarize in three simple words the daily perception of the users of **RESISTA** branded products.

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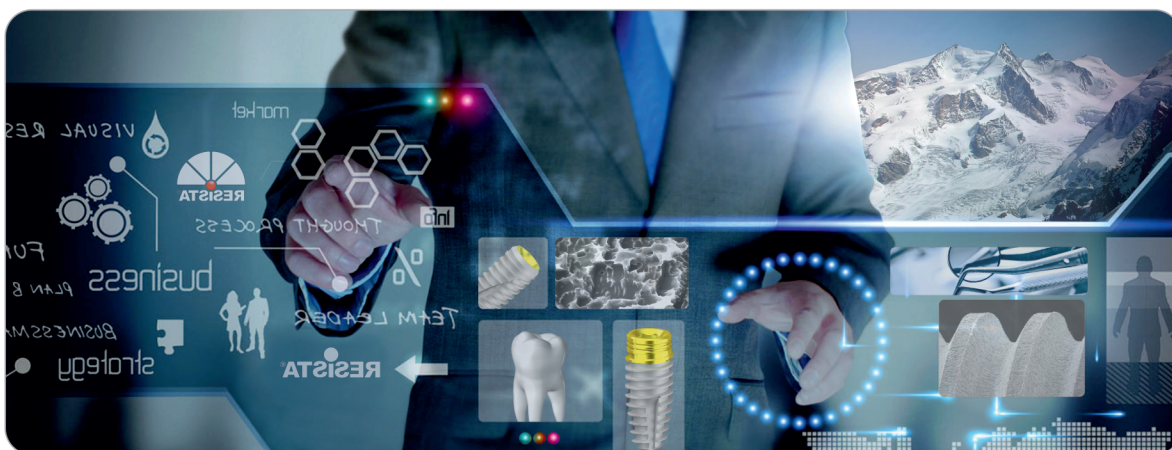
RESEARCH & DEVELOPMENT - The R&D Resista Group's program stems from indications from the clinical world combined with our experience in the field of implantable medical devices.

Due to the use of 3D modeling programs and sophisticated computer systems, we simulate the final geometries and the design, developing the rapid prototyping phases of the projects.

INNOVATIONS - Resista Group is structured to meet all the needs of the metal-implant with an advanced conception technology, prototyping and final product engineering to arrive at innovative solutions. Validation is carried out in collaboration with accredited Research Institutes following the International Standards Requirements.

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PROJECT AND PROCESS QUALITY - ICIM Spa certified **Ing. C. A. Issoglio & C. S.r.l.** in accordance with **UNI EN ISO 9001** and **UNI EN ISO 13485** in compliance with all current regulations relating to the products and services offered. Every employee covering different roles (technical, engineering, administrative, commercial, mechanic operators, etc.), follows guidelines and objectives for a single purpose: continuous product improvement.

MAIN OBJECTIVE - Customer satisfaction is our objective. The winning factor is the company's ability to quickly resolve any type of request. The technical support provided by Product Specialists, availability, efficiency and friendliness are the strengths of our structure.

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BIOMECHANICAL CONCEPTS - The **Resista Group** method success is the result of several product design studies and the high technology applied in the production. The precision and reliability of the implant has a strong influence on the long term clinical success.


The manufacturing process is carried out by high qualified operators, skilled enough to obtain optimal and reproducible results with the use of electronically controlled CNC machinery. The machinery tolerances, especially in **prosthetic connections**, are our **gold standard** and they are made between **7 - 10 µm**.

MATERIALS - The implants are made in **Titanium ASTM Gr.4 Cold Worked** (ISO 5832/2), the prosthetic screws and the prosthetic components are made in **Titanium ASTM Gr.5 ELI** (ISO 5832/2).

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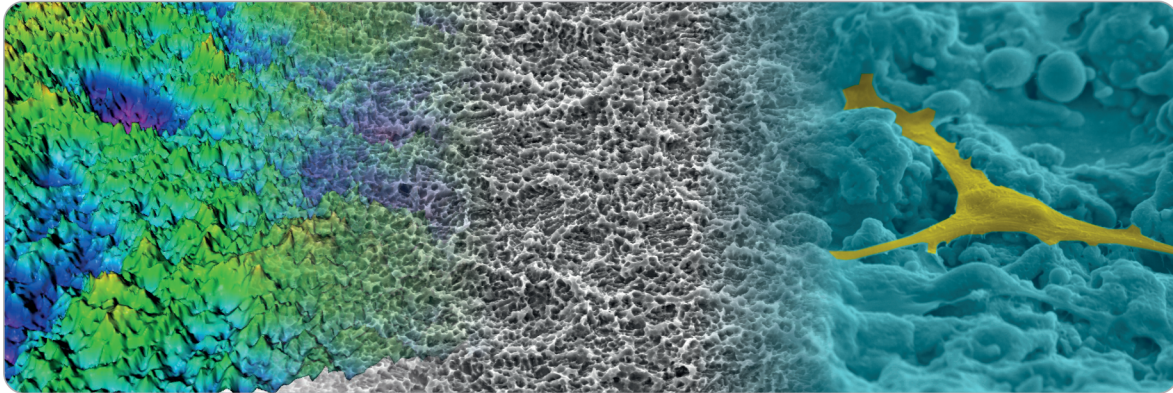
SURFACE TREATMENTS - The non-contaminant micro-subtraction process, **DAE** (Double Acid Etching), changes the implant micro roughness (**Ra, Rq**) and the surface texture maximizing the wettability and biomimetic properties. The new surface treatment and decontamination process are patented by **Nobil Bio Ricerche** improving the bio-chemical properties of the implants.

QUALITY SYSTEM (QS) - UNI EN ISO 9001 and **UNI EN ISO 13485** in compliance of Medical Devices Directives. We are qualified in the design and in the production management of dental implants, dental prosthesis, intraligamental anesthesia syringes and abrasive discs.

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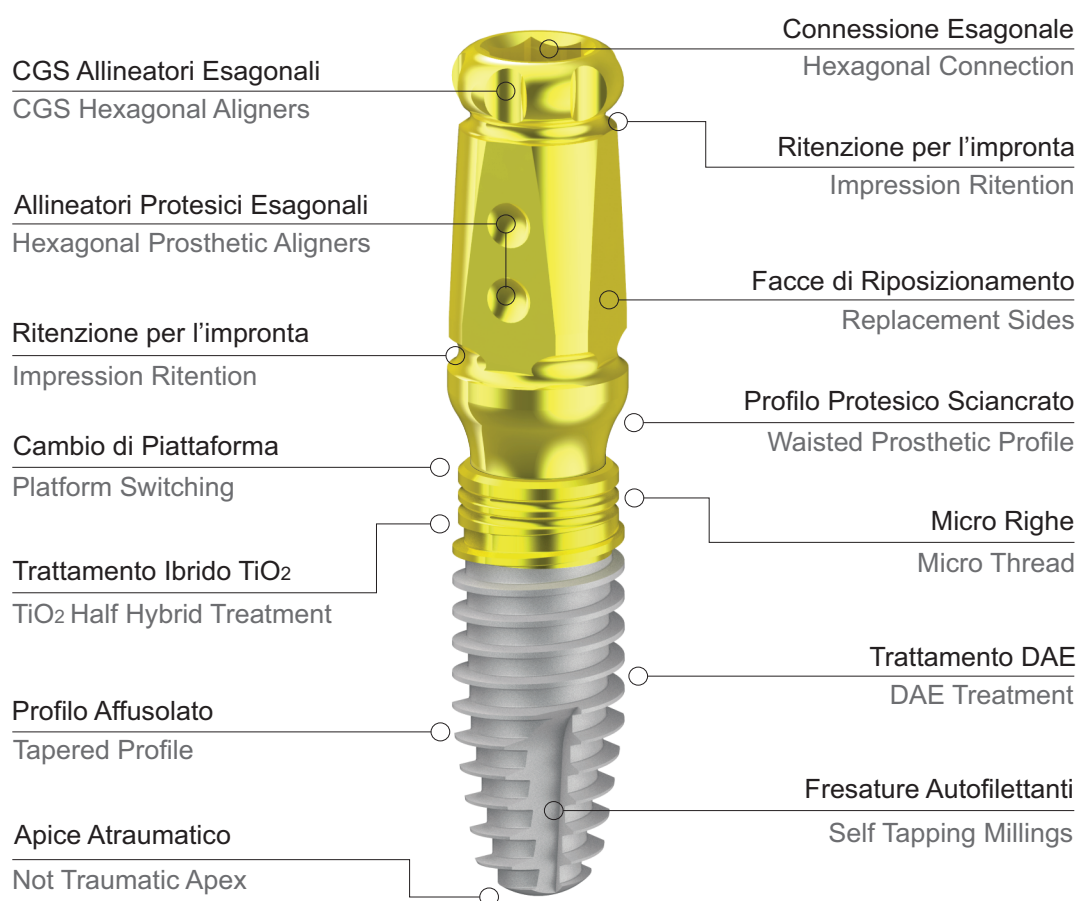
TITANIUM GR 5 PRE-ASSEMBLED MOUNTER - One of the strong points of our **Resista** line is the Multifunction Mounter Device that is used in the **4** following different applications:

1 Implant pick up and screwing - **2** Closed Tray Impression Transfer - **3** Definitive Abutment up to 13° angular correction - **4** Mounter for Computer Guided Surgery

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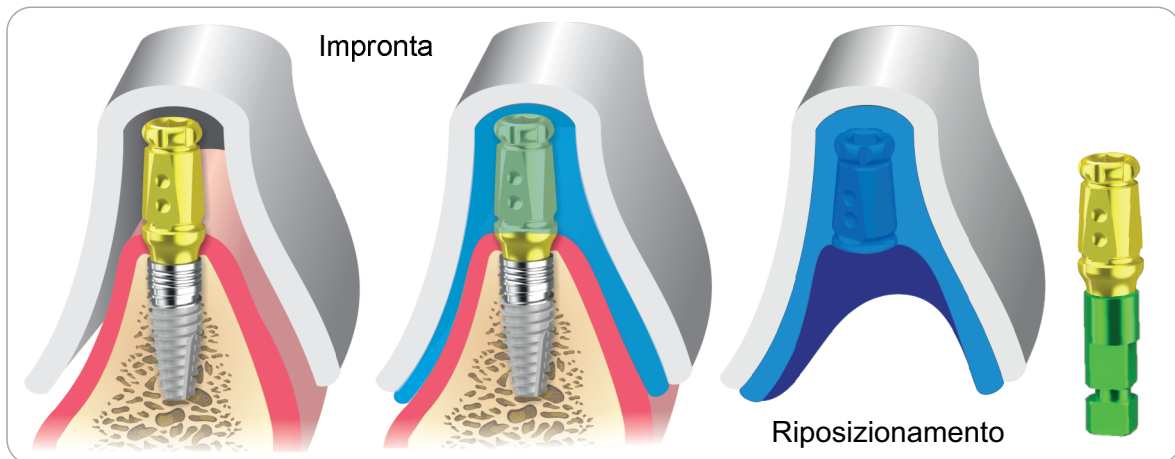TRANSFER MOUNTER - Transfer for **Closed Tray Technique** (Elastic tear impression technique).

It is a simple, fast and effective technique, suitable for all types of prosthetic projects. The Transfer Mounter can be repositioned in the impression only in one position, defined by the different sides and retentions arranged on the surface.


ABUTMENT MOUNTER - Definitive Prosthetic Abutment complete with definitive through screw.

The Gr 5 Titanium Abutment Mounter is certified to be used as a permanent abutment, can be milled in height and angle up to 13°, using Titanium laboratory burs.

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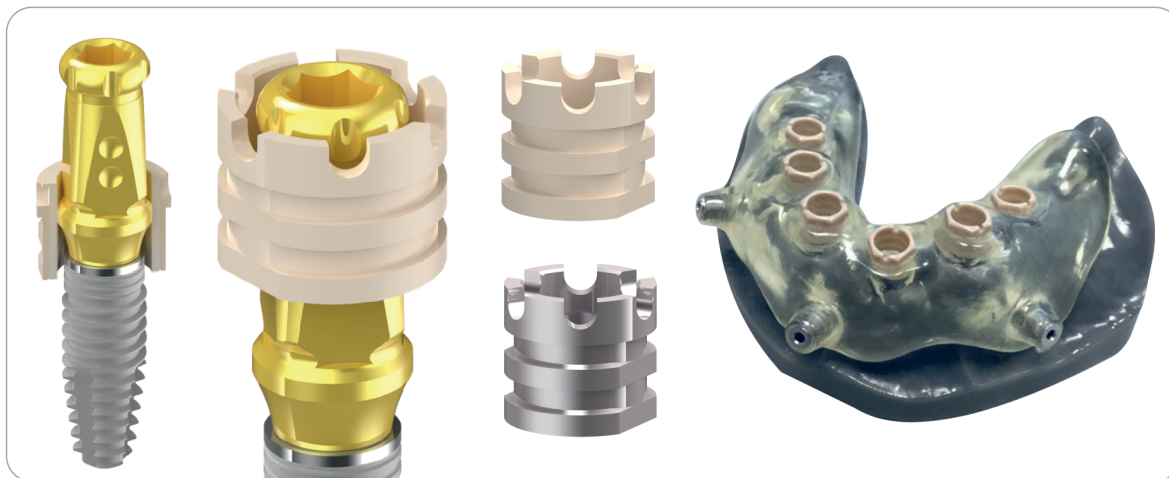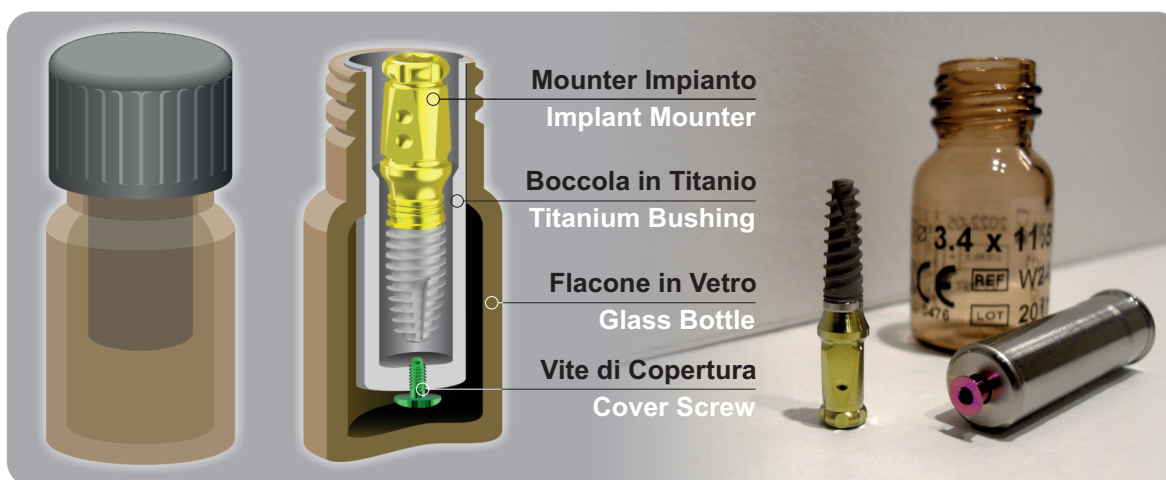
CGS MOUNTER - Mount suitable for implants placement by **Computer Guided Surgery**.

The cylindrical diameter of the Mounter glides on the axis inside the Bushings for Computer Guided Surgery, whereas the 6 hexagonal alignment reference points (notches) define the final position.

THE PACKAGING - The **Titanium Bushing** contains the Mounter and the implant and guarantees the quality of the surface during the removal surgical step, avoiding the contamination due to the scratching. The mounter and implant are contained within a **Gr. 4 Titanium Bushing**, which guarantees surface quality over time and during surgical removal, preventing contamination.

The **Glass Bottle** is the most suitable material for preserving the surface chemical properties of the titanium implant. This prevents the implant from absorbing the shielding hydrocarbon atoms, released by the plastic packaging when irradiated with Gamma or Beta rays.

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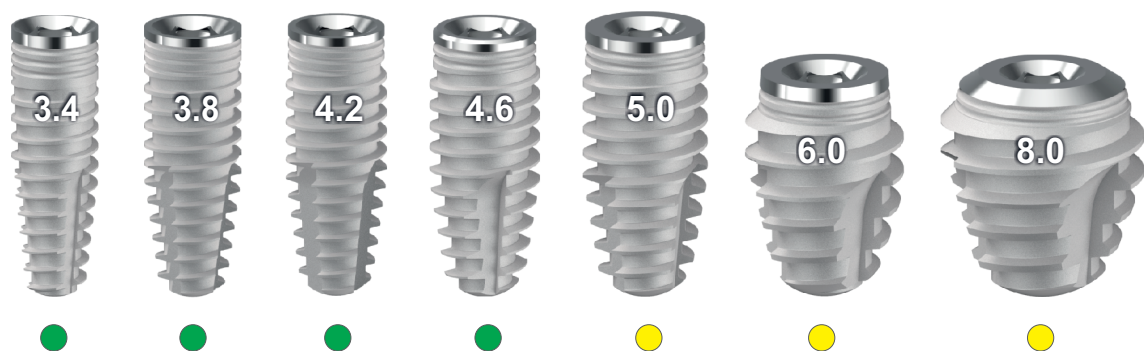
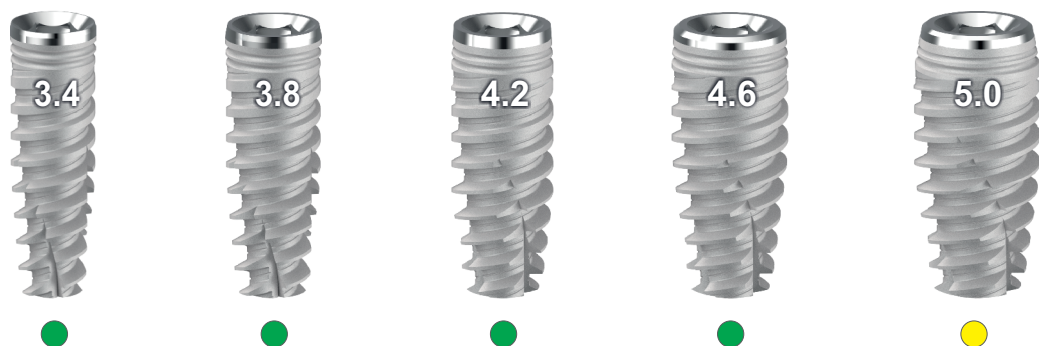
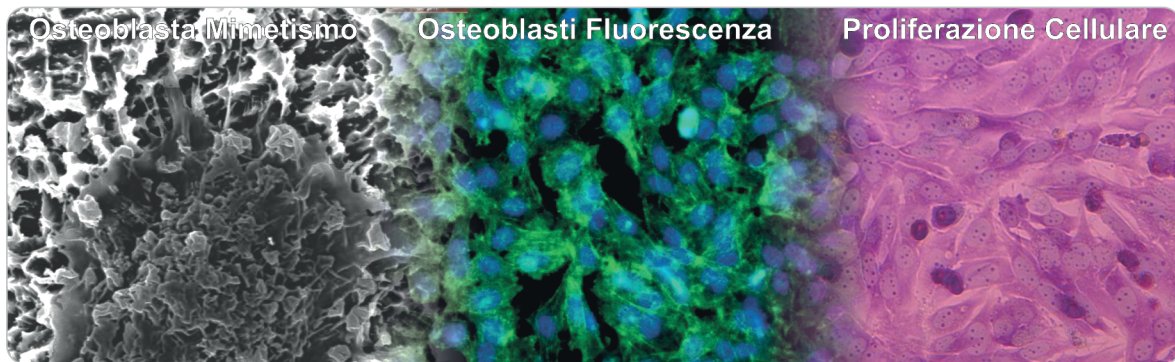
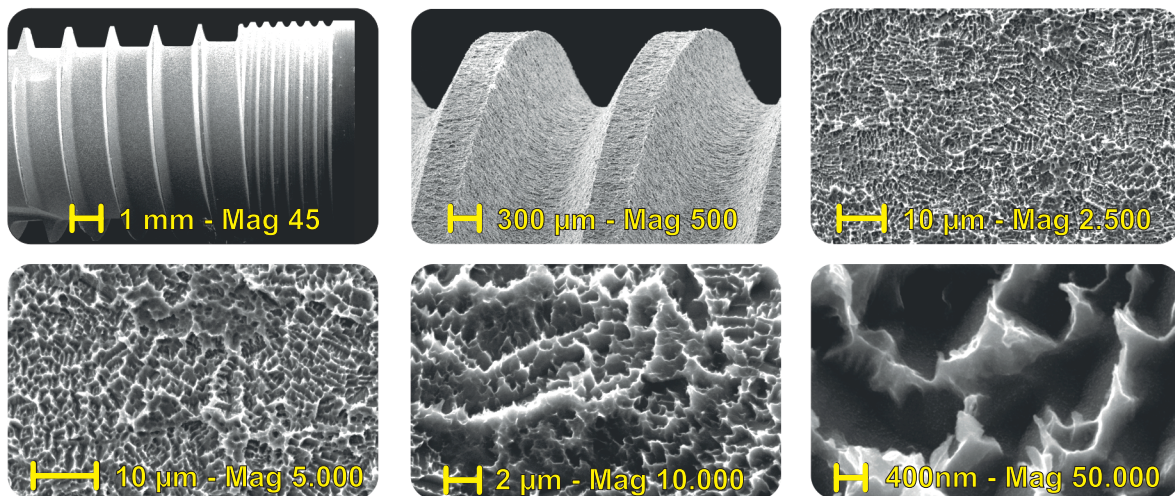
The Micro-Nano Rough DAE surface treatment speed up bone healing processes

1. Removes the manufacturing organic residuals along with Argon Plasma Cleaning
2. Smooths the edge and eliminates the micro-defects
3. Increases surface and wettability (hydrophilicity surface), improving the first fibrin bridges adhesion
4. Increases protein adhesion
5. Maximizes cellular adhesion through the Micro-Nano roughness suitable for actin filaments (phyllopods) anchoring
6. Changes the Titanium surface chemical characteristics improving the biocompatibility and increasing the cellular proliferation and vitality

Surface topography and chemistry interact with the cellular differentiation processes


The implant surface uses the microtopography as a communication language with the host tissue cells. The Plasma Argon Cold Cleaning with a controlled-environment packaging at a strict protocol procedures, play a key-role in the control of adhered endotoxins (main immunological response player to implant surfaces).

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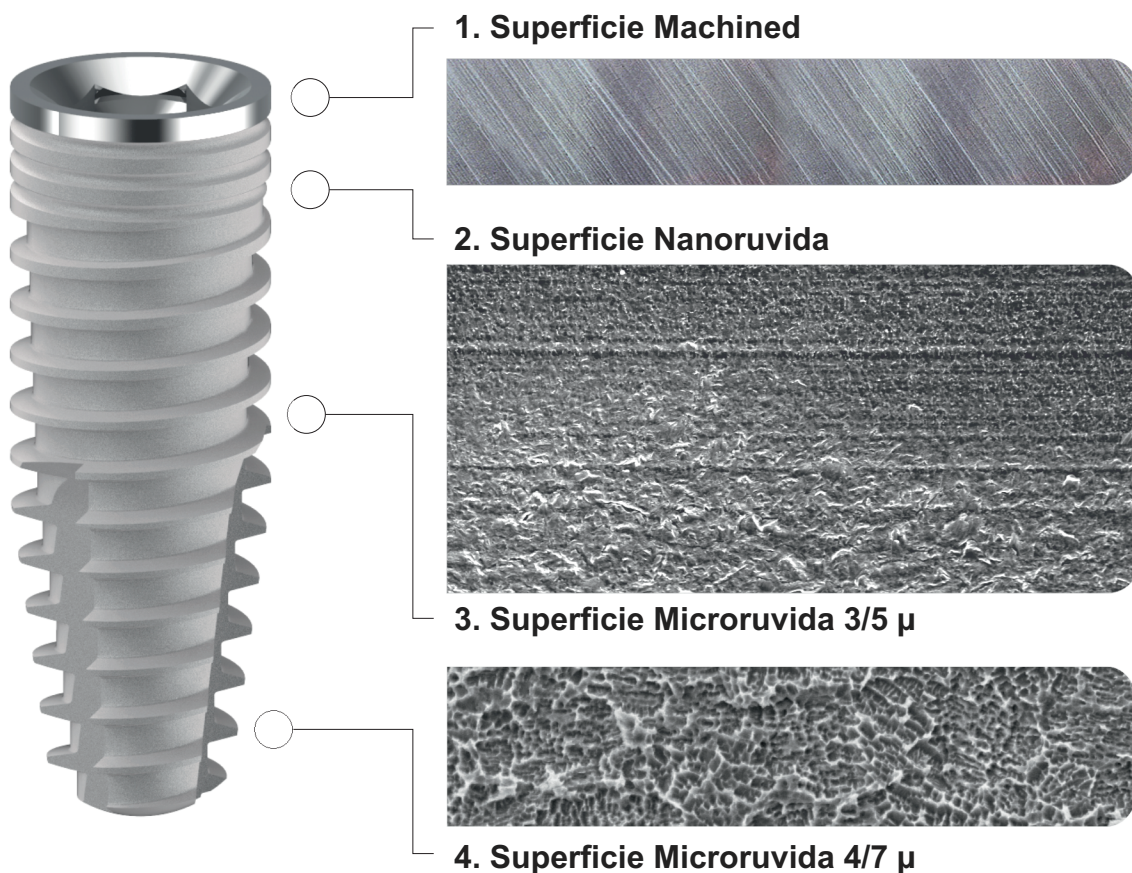VEGF/VEGF-R/RUNX2 Upregulation in Human Periodontal Ligament Stem Cells Seeded on Dual Acid Etched Titanium Disk

Francesca Diomedè, Guya Diletta Marconi, Marcos F. X. B. Cavalcanti, Jacopo Pizzicannella, Sante Donato Pierdomenico, Luigia Fonticoli, Adriano Piattelli and Oriana Trubiani

Materials 2020, 13, 706; doi:10.3390/ma13030706


SURFACE TREATMENT DISTRIBUTION ON IMPLANT'S BODY

All the **Resista**'s implants have a surface treatment that is locally differentiated. The logic of surface treatment distribution is dictated by the endo-osseous implant position which requires different characteristics dependin on the muco-gingival component presence around.

1. **Machined Surface:** Brushed, rectify and protected during the roughening process
2. **Nano-rough Surface:** Speedy treatment, nano-rough with average Ra less than one micron
3. **Micro-rough Surface:** 3/5 μ : DAE Treatment, micro-rough with average Ra between 3/5 μ
4. **Micro-rough Surface:** 4/7 μ : DAE Slow Treatment, micro-rough with average Ra between 4/7 μ

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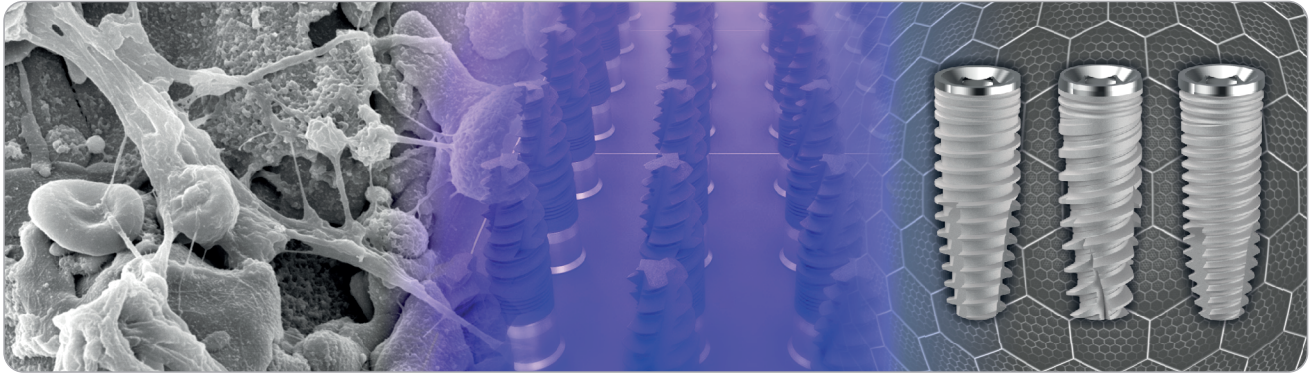The Bacterial Anti-Adhesive Activity of Double-Etched Titanium (DAE) as a Dental Implant Surface

Morena Petrini, Alessandra Giuliani, Emanuela Di Campi, Silvia Di Lodovico, Giovanna Iezzi, Adriano Piattelli and Simonetta D'Ercole¹
[International Journal of Molecular Sciences 2020, 21, 8315; doi:10.3390/ijms21218315](#)

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THE RESISTA'S IMPLANT LINE STRONG POINTS

- 1. Macro Geometry:** A vast range of implant shapes for different surgical needs.
- 2. Micro Geometry:** Last generation differentiated surface treated (DAE Micro and Nano rough), Certified Decontamination in Argon Cold Plasma Reactor with white clean room packaging.
- 3. Multifunctional Mounter:** 4 components in a single article (Mounter, Transfer, Abutment and Mounter for Computer Guided Surgery) to maximize result and minimize costs.
- 4. Prosthetic Parts:** A very wide range of prosthetic variants with 7 micron machining tolerances on the hexagon connection, with emerging curvilinear profiles, flat and conical connections, platform switching, golden titanium screws and PVD TiN coatings (more aesthetic and less retention of dental plaque).
- 5. Surgical Instruments:** Surgical Drills **3-Tech**, maximum efficiency, minimum trauma, smaller friction coefficient, **PRO MSD** technology for counterclockwise bone densifying application long, lasting resistance, perfect visibility of laser marking and surgical millimeters stop.

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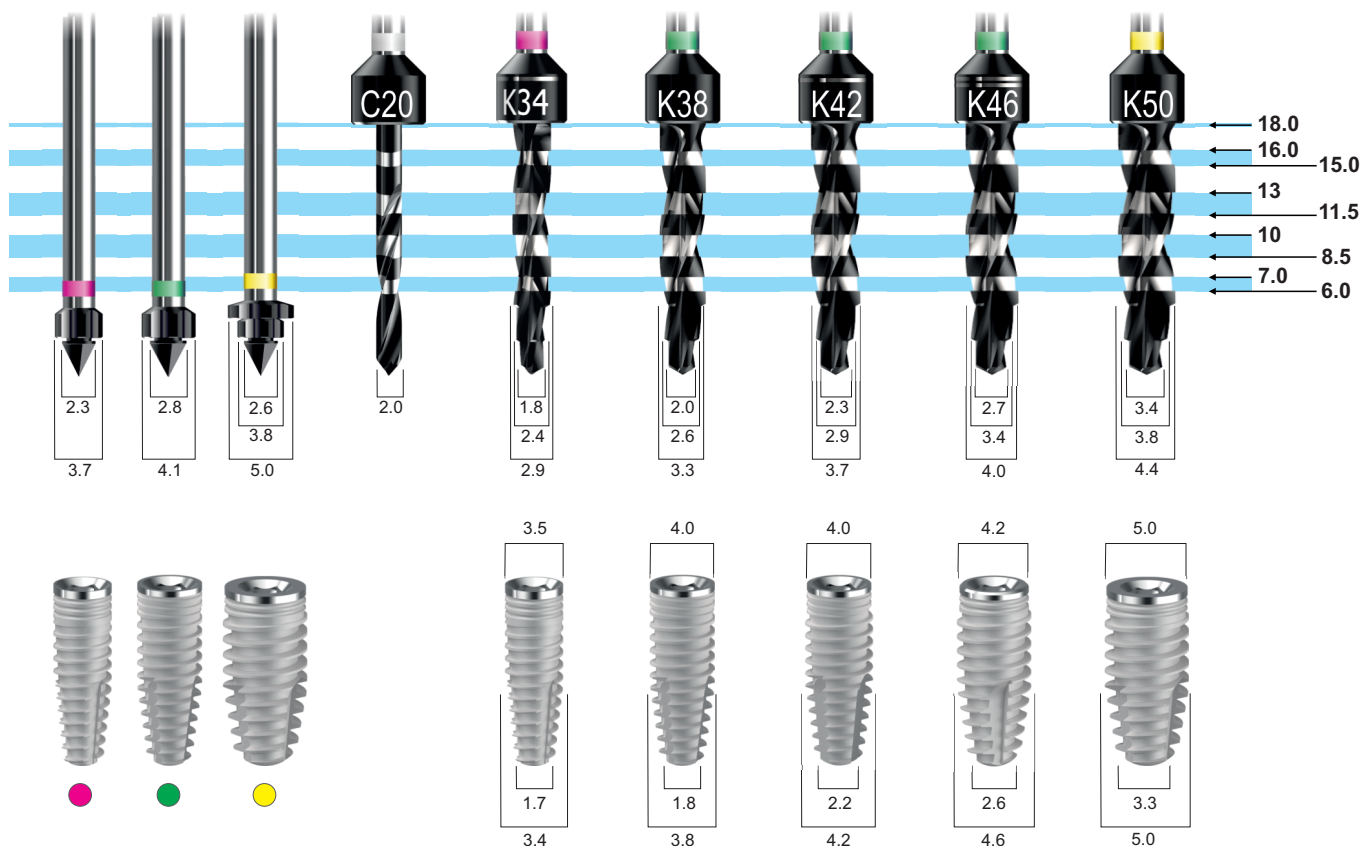
✓ Various bio-mechanical factors affecting heat generation during osteotomy preparation: A systematic review

Chirag J Chauhan¹, Darshana N Shah¹, Foram B Sutaria¹

Indian J Dent Res. Jan-Feb 2018;29(1):81-92. doi: 10.4103/ijdr.IJDR_729_16.

LOGIC RELATION BETWEEN DRILL AND IMPLANT CORE

The rotary instruments design and shape is based on the logical coherence between the implants size, in the external dimensions, internal core and their penetration properties.

The **IC line Cylindrical Implants** have a tapered apical portion 3mm long.

The **IK line Conical Implants** have a tapered apical portion 5mm long.

The drills have respectively 2/3 apical sections, suitable for the implant tunnel preparation, with respect to the implant geometry, maintaining the best bone implant contact (BIC).

Each step drill is equipped with a sharp cutting edge to facilitate penetration, minimizing friction and overheating.

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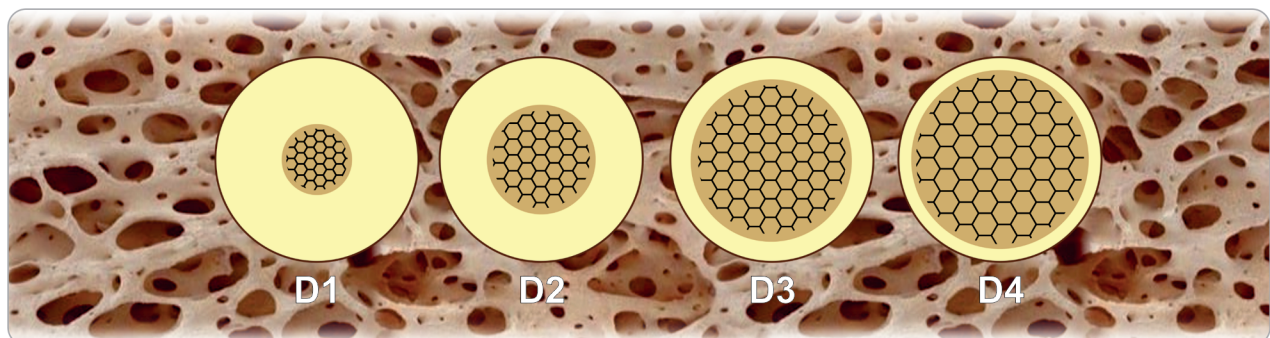
✓ **Bone classification: clinical-histomorphometric comparison**

Trisi P, Rao W. - Clin Oral Implants Res. 1999 Feb;10(1):1-7. doi: 10.1034/j.1600-0501.1999.100101.x.

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THE CONDITIONING OF SURGICAL PROTOCOLS

The differentiated implant geometries in the macro/micro shape provide for different applications, both according to the surgical/prosthetic requirements (immediate loading, early loading or postponed loading) and according to the conditions of the bone/anatomical component of the patient who needs rehabilitation. However, surgical maneuvers and alternative protocols are often required to improve implant performance, according to the morphological characteristics of the residual bone structures.

The success of these maneuvers, that is often the result of a customized tools and protocols combination, represents the different response to the problem in relation to the personal operating skills and the presence of tools suitable for their application.

MECHANICAL OSSEODENSIFICATION IN ANTI-CLOCKWISE ROTATION

The osseodensification using rotary instruments is a recent surgical technique for the implant site preparation that can be associated with different protocols; these protocols can be used in those particular anatomical conditions such as poor bone quality and insufficient vertical / horizontal dimensions.

This non-subtractive osteotomic bone-condensation approach produces an increase in peri-osteotomic bone density, saving the tissue and increasing the primary implant stability.

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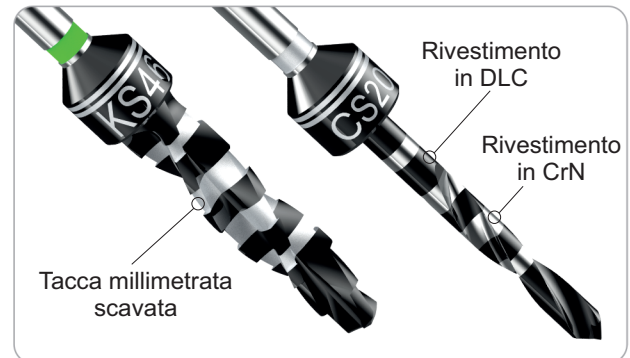
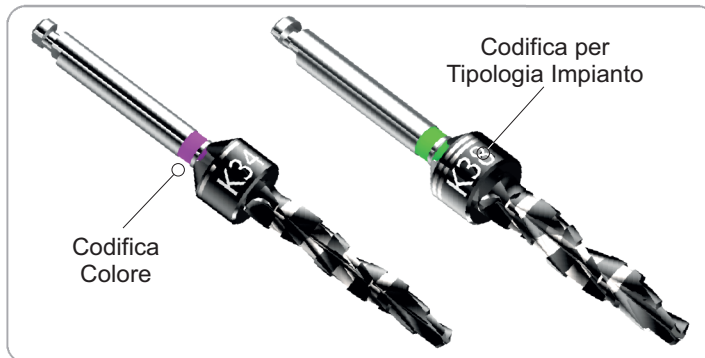
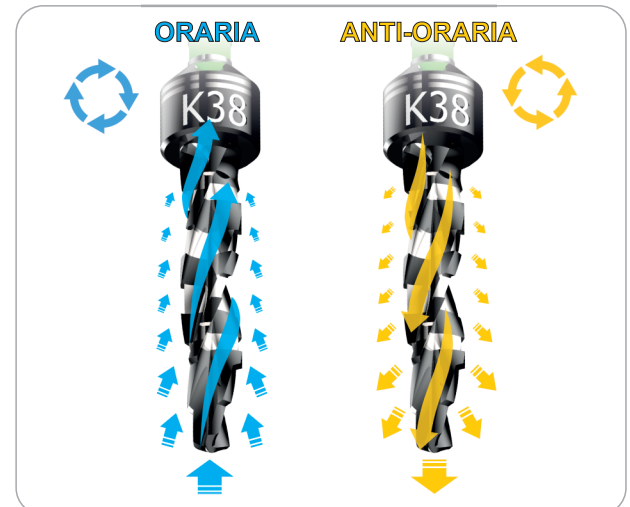
RESISTA PRO MSD SURGICAL DRILLS - Modular Surgical Drilling

The Modular Surgical Drills with variable geometry are 3-section helicoidal drills with progressive radius rake, that can be used clockwise and anticlockwise.

The diameter of each drill changes according to a constant ratio of 0.4mm (30/34/38/42/46/50): this allows the operator to choose the drill according to the bone quality (over-preparation or under-preparation). All the drills are coated with a DLC (Diamond Like Carbon) synthetic diamond solid lubricant that maximizes performance in terms of mechanical strength and friction reduction.


The **anticlockwise rotation**, reversing the involved forces, generates 3 different effects on the implant tunnel perforation, which can revolutionize the known milling logic.

- 1) Anterior and lateral thrust of the bone removed by the tip, and of the liquids as well, blood and physiological water.
- 2) Ejection and "hammer effect" of the drill that produce an improvement in vertical control.
- 3) Reduction of cutting efficiency to protect sensitive anatomical parts.

INTERNAL HEXAGON UNIVERSAL CONNECTION

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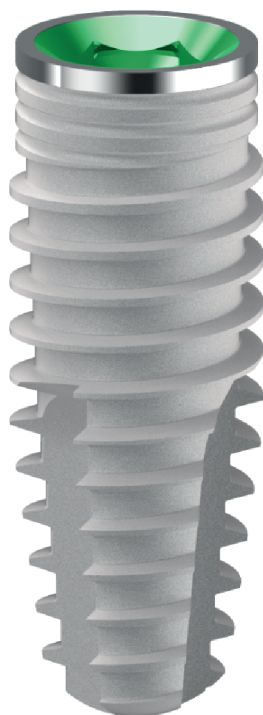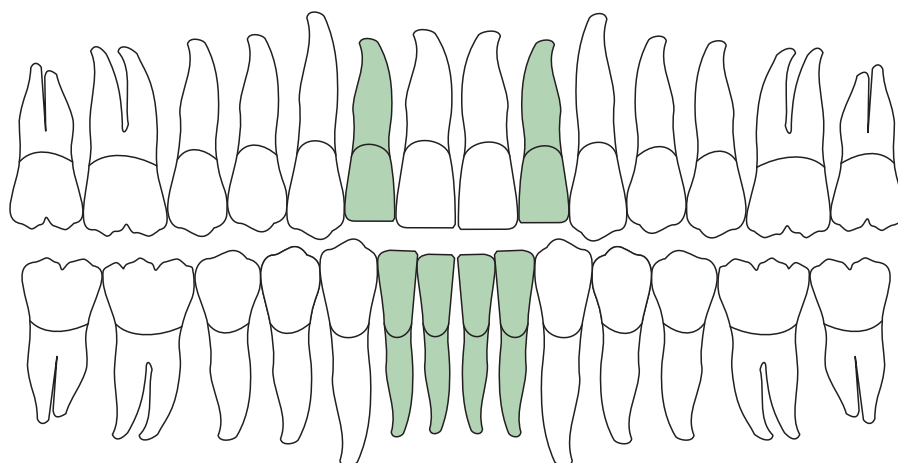
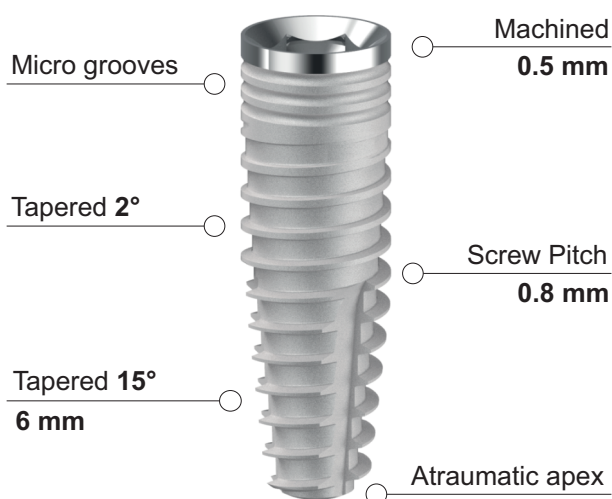
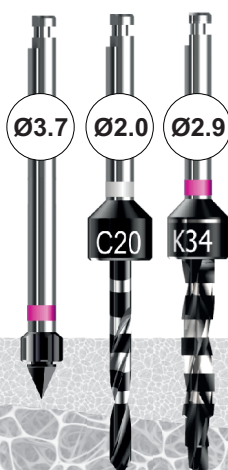
Recommended milling burs rotation speed 500 - 800 rpm

3.4 Full Treatment - With Mounter

CODE	IMPLANT	LENGTH	PLATFORM	APEX
IK 3407	Ø 3.4	7.0 mm	Ø 3.5 CPS	Ø 1.8
IK 3408	Ø 3.4	8.5 mm	Ø 3.5 CPS	Ø 1.8
IK 3410	Ø 3.4	10 mm	Ø 3.5 CPS	Ø 1.8
IK 3411	Ø 3.4	11.5 mm	Ø 3.5 CPS	Ø 1.8
IK 3413	Ø 3.4	13 mm	Ø 3.5 CPS	Ø 1.8
IK 3415	Ø 3.4	15 mm	Ø 3.5 CPS	Ø 1.8

3.4 Half Treatment **Anatase TiO₂** - With Mounter

CODE	IMPLANT	LENGTH	PLATFORM	APEX
IK 3408 HT	Ø 3.4	8.5 mm	Ø 3.5 CPS	Ø 1.8
IK 3410 HT	Ø 3.4	10 mm	Ø 3.5 CPS	Ø 1.8
IK 3411 HT	Ø 3.4	11.5 mm	Ø 3.5 CPS	Ø 1.8
IK 3413 HT	Ø 3.4	13 mm	Ø 3.5 CPS	Ø 1.8
IK 3415 HT	Ø 3.4	15 mm	Ø 3.5 CPS	Ø 1.8

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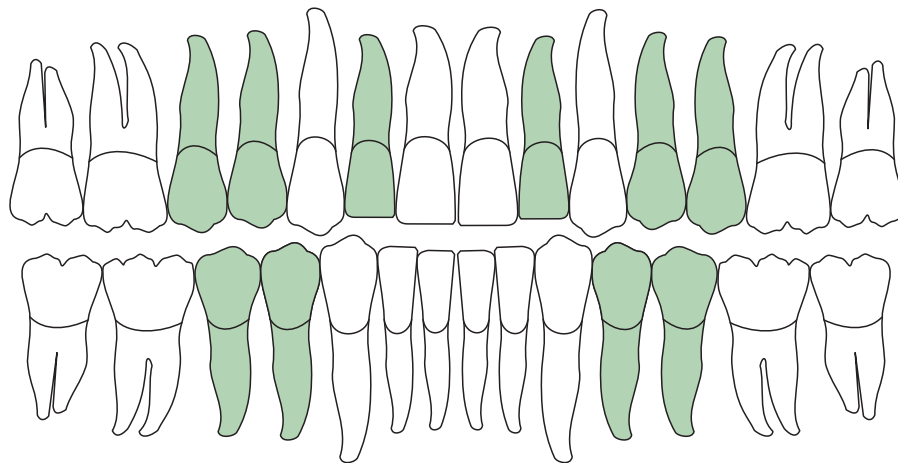
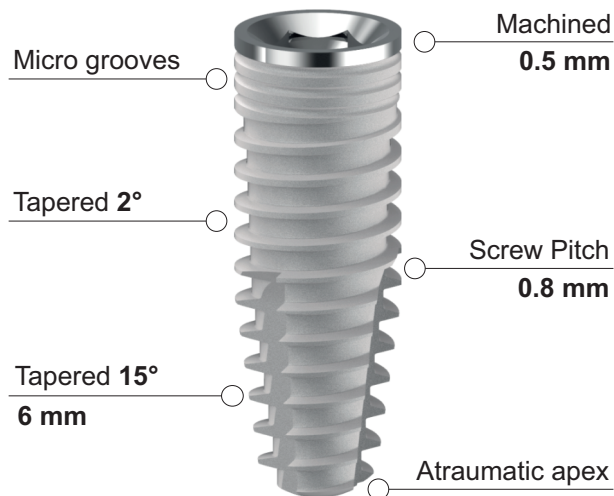
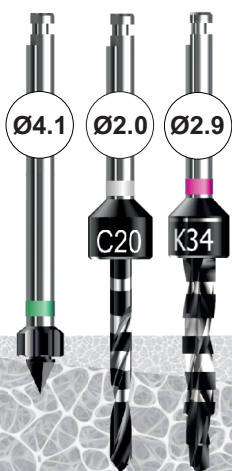
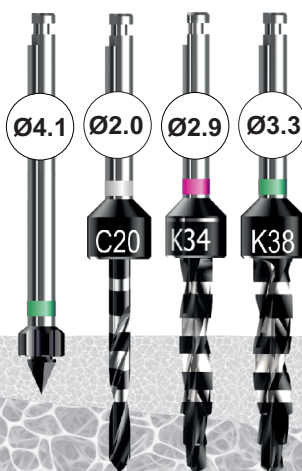
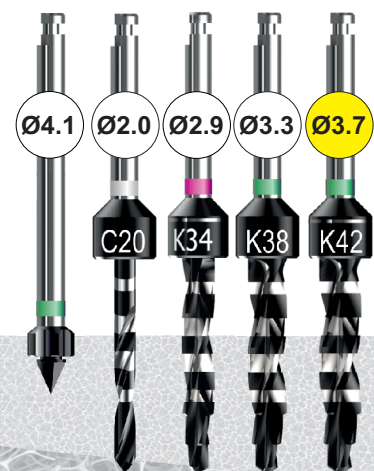
Recommended milling burs rotation speed 500 - 800 rpm

3.8 Full Treatment - With Mounter

CODE	IMPLANT	LENGTH	PLATFORM	APEX
IK 3808	Ø 3.8	8.5 mm	Ø 4.0	Ø 1.9
IK 3810	Ø 3.8	10 mm	Ø 4.0	Ø 1.9
IK 3811	Ø 3.8	11.5 mm	Ø 4.0	Ø 1.9
IK 3813	Ø 3.8	13 mm	Ø 4.0	Ø 1.9
IK 3815	Ø 3.8	15 mm	Ø 4.0	Ø 1.9
IK 3818	Ø 3.8	18 mm	Ø 4.0	Ø 1.9

3.8 Half Treatment **Anatase TiO₂** - With Mounter

CODE	IMPLANT	LENGTH	PLATFORM	APEX
IK 3808 HT	Ø 3.8	8.5 mm	Ø 4.0	Ø 1.9
IK 3810 HT	Ø 3.8	10 mm	Ø 4.0	Ø 1.9
IK 3811 HT	Ø 3.8	11.5 mm	Ø 4.0	Ø 1.9
IK 3813 HT	Ø 3.8	13 mm	Ø 4.0	Ø 1.9
IK 3815 HT	Ø 3.8	15 mm	Ø 4.0	Ø 1.9

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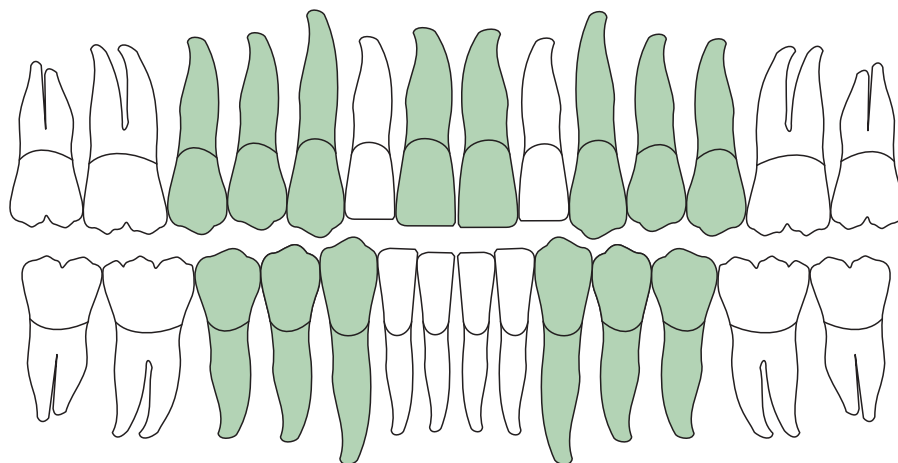
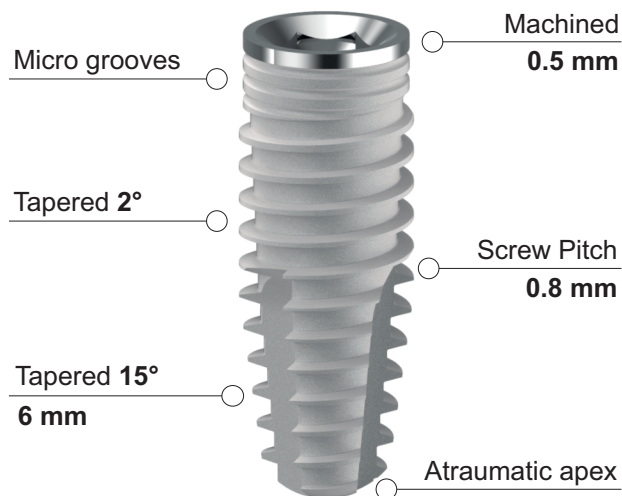
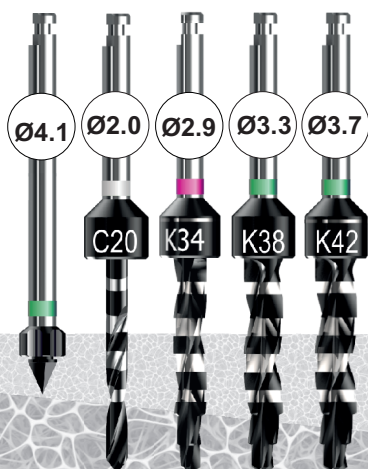
Recommended milling burs rotation speed 500 - 800 rpm

4.2 Full Treatment - With Mounter

CODE	IMPLANT	LENGTH	PLATFORM	APEX
IK 4207	Ø 4.2	7.0 mm	Ø 4.0	Ø 2.3
IK 4208	Ø 4.2	8.5 mm	Ø 4.0	Ø 2.3
IK 4210	Ø 4.2	10 mm	Ø 4.0	Ø 2.3
IK 4211	Ø 4.2	11.5 mm	Ø 4.0	Ø 2.3
IK 4213	Ø 4.2	13 mm	Ø 4.0	Ø 2.3
IK 4215	Ø 4.2	15 mm	Ø 4.0	Ø 2.3

4.2 Half Treatment **Anatase TiO₂** - With Mounter

CODE	IMPLANT	LENGTH	PLATFORM	APEX
IK 4208 HT	Ø 4.2	8.5 mm	Ø 4.0	Ø 2.3
IK 4210 HT	Ø 4.2	10 mm	Ø 4.0	Ø 2.3
IK 4211 HT	Ø 4.2	11.5 mm	Ø 4.0	Ø 2.3
IK 4213 HT	Ø 4.2	13 mm	Ø 4.0	Ø 2.3
IK 4215 HT	Ø 4.2	15 mm	Ø 4.0	Ø 2.3

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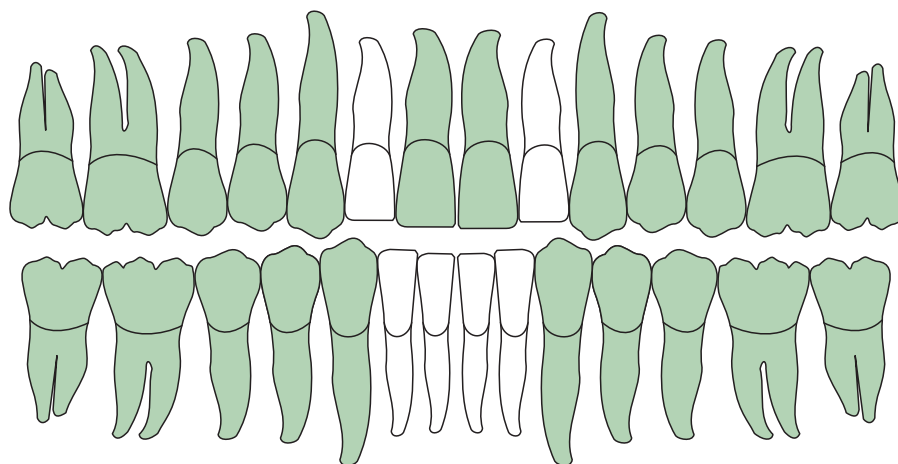
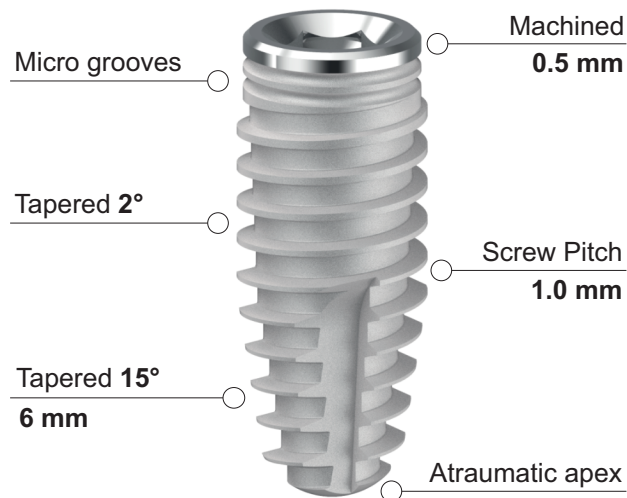
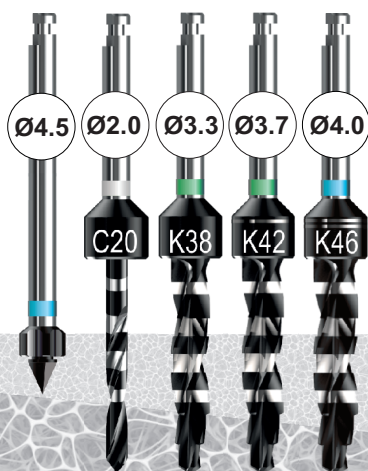
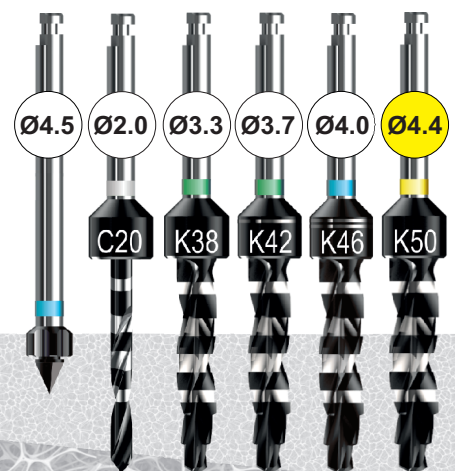
Recommended milling burs rotation speed 500 - 800 rpm

4.6 Full Treatment - With Mounter

CODE	IMPLANT	LENGTH	PLATFORM	APEX
IK 4607	Ø 4.6	7.0 mm	Ø 4.2	Ø 2.5
IK 4608	Ø 4.6	8.5 mm	Ø 4.2	Ø 2.5
IK 4610	Ø 4.6	10 mm	Ø 4.2	Ø 2.5
IK 4611	Ø 4.6	11.5 mm	Ø 4.2	Ø 2.5
IK 4613	Ø 4.6	13 mm	Ø 4.2	Ø 2.5
IK 4615	Ø 4.6	15 mm	Ø 4.2	Ø 2.5

4.6 Half Treatment **Anatase TiO₂** - With Mounter

CODE	IMPLANT	LENGTH	PLATFORM	APEX
IK 4608 HT	Ø 4.6	8.5 mm	Ø 4.2	Ø 2.5
IK 4610 HT	Ø 4.6	10 mm	Ø 4.2	Ø 2.5
IK 4611 HT	Ø 4.6	11.5 mm	Ø 4.2	Ø 2.5
IK 4613 HT	Ø 4.6	13 mm	Ø 4.2	Ø 2.5
IK 4615 HT	Ø 4.6	15 mm	Ø 4.2	Ø 2.5

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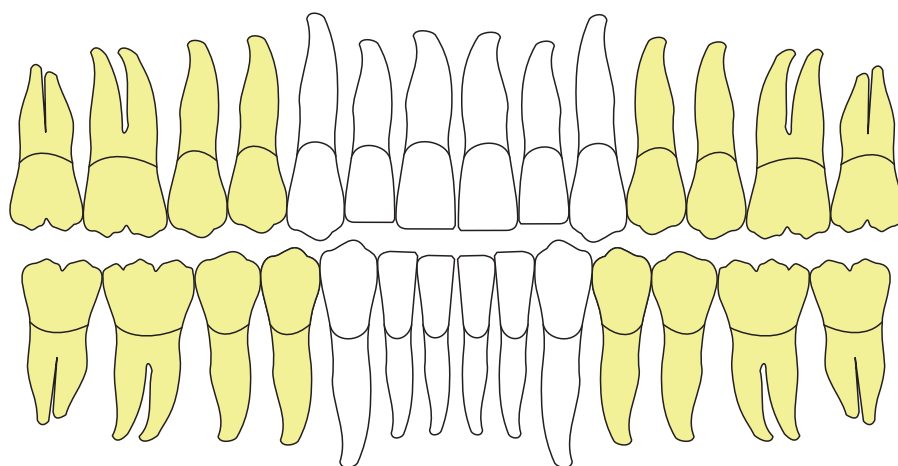
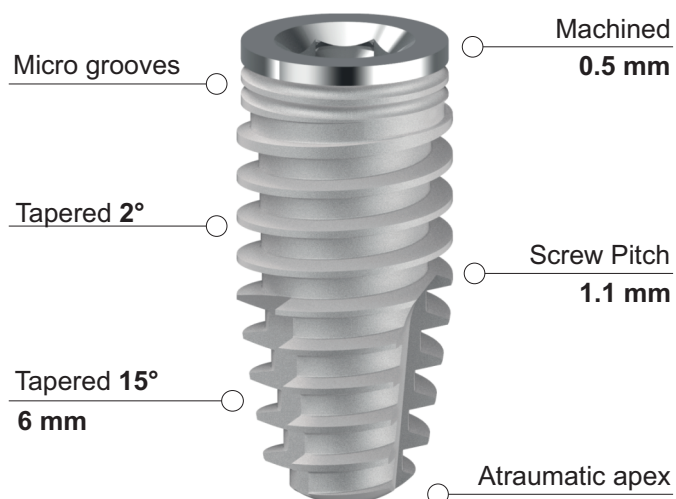
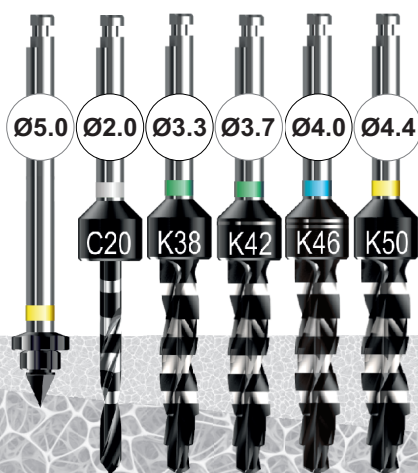
Recommended milling burs rotation speed 500 - 800 rpm

5.0 Full Treatment - With Mounter

CODE	IMPLANT	LENGTH	PLATFORM	APEX
IK 5007	Ø 5.0	7.0 mm	Ø 5.0	Ø 2.7
IK 5008	Ø 5.0	8.5 mm	Ø 5.0	Ø 2.7
IK 5010	Ø 5.0	10 mm	Ø 5.0	Ø 2.7
IK 5011	Ø 5.0	11.5 mm	Ø 5.0	Ø 2.7
IK 5013	Ø 5.0	13 mm	Ø 5.0	Ø 2.7

5.0 Half Treatment **Anatase TiO₂** - With Mounter

CODE	IMPLANT	LENGTH	PLATFORM	APEX
IK 5008 HT	Ø 5.0	8.5 mm	Ø 5.0	Ø 2.7
IK 5010 HT	Ø 5.0	10 mm	Ø 5.0	Ø 2.7
IK 5011 HT	Ø 5.0	11.5 mm	Ø 5.0	Ø 2.7
IK 5013 HT	Ø 5.0	13 mm	Ø 5.0	Ø 2.7


5.0 Full Treatment - Mounterless


CODE	IMPLANT	LENGTH	PLATFORM	APEX
IK 5007 ML	Ø 5.0	7.0 mm	Ø 5.0	Ø 2.7
IK 5008 ML	Ø 5.0	8.5 mm	Ø 5.0	Ø 2.7
IK 5010 ML	Ø 5.0	10 mm	Ø 5.0	Ø 2.7
IK 5011 ML	Ø 5.0	11.5 mm	Ø 5.0	Ø 2.7
IK 5013 ML	Ø 5.0	13 mm	Ø 5.0	Ø 2.7

IK 60 - IK 80

X-LARGE TAPERED SHAPE Internal Hexagon

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
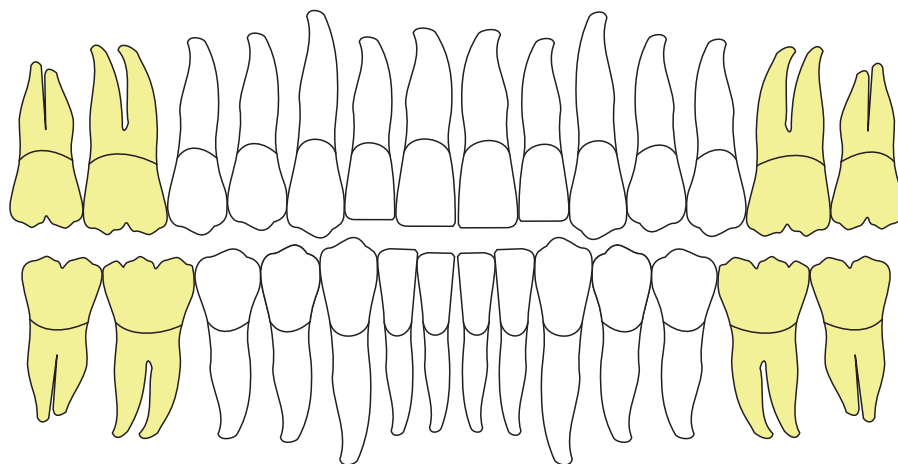
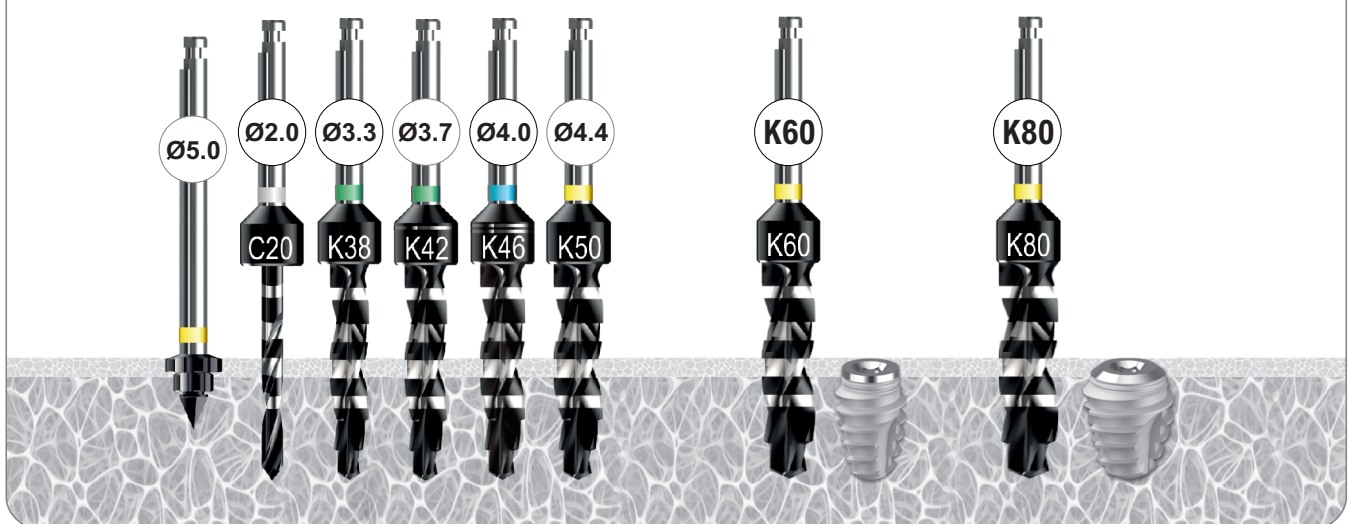
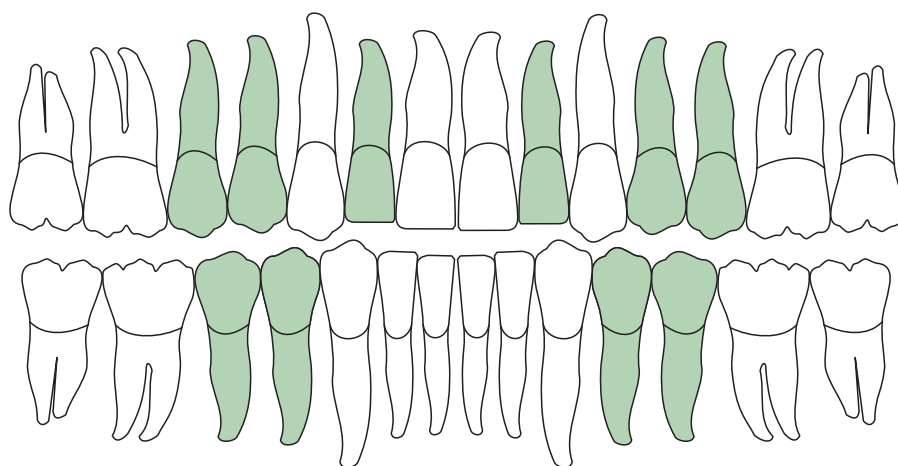
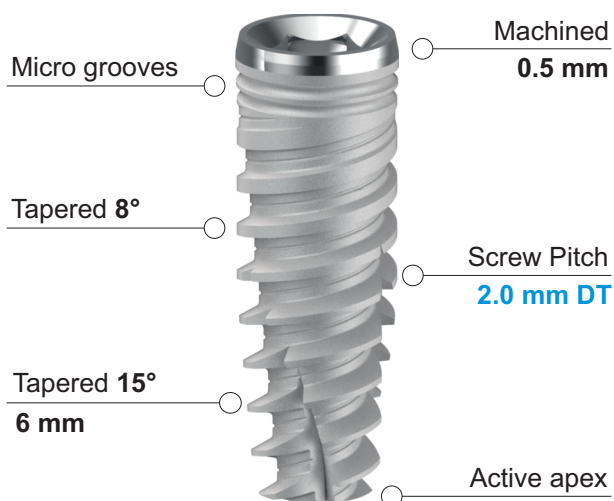
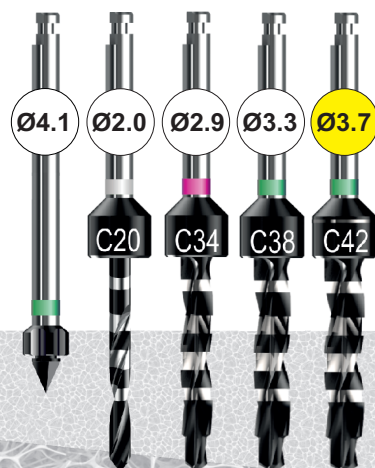
THE RIGHT POST-EXTRACTION SURGICAL TECHNIQUE

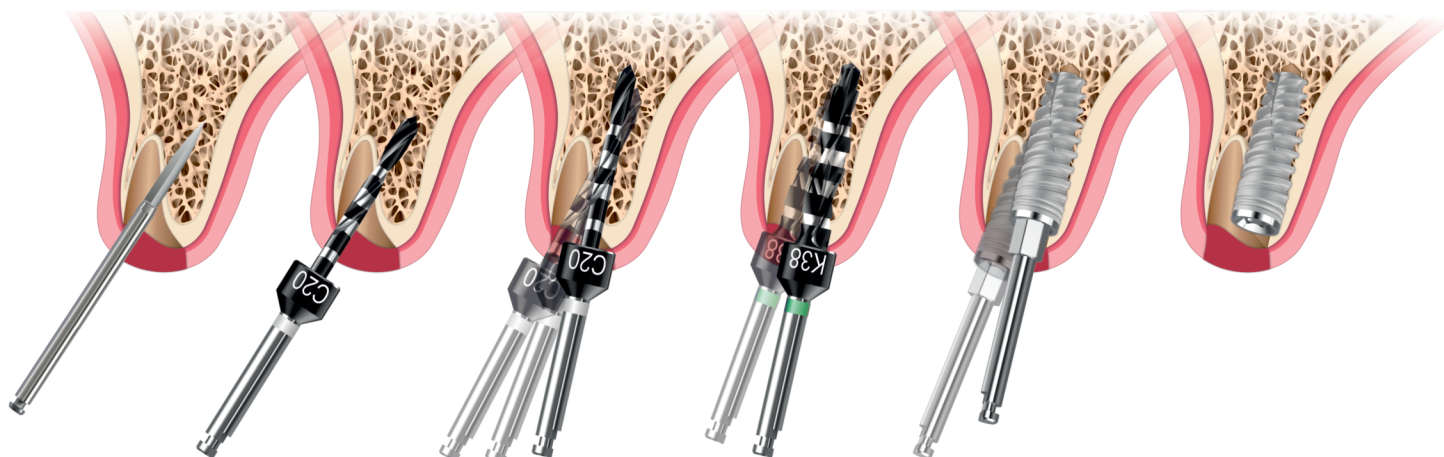IA 46


ACTIVE BLADE SHAPE Internal Hexagon

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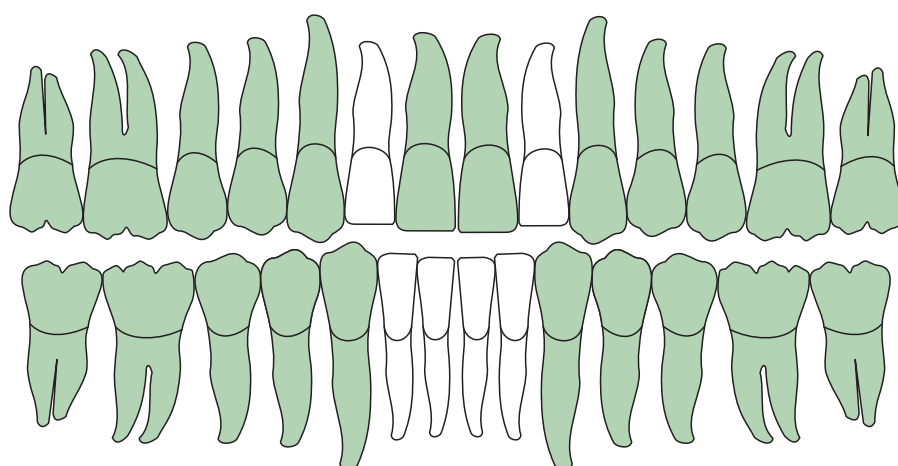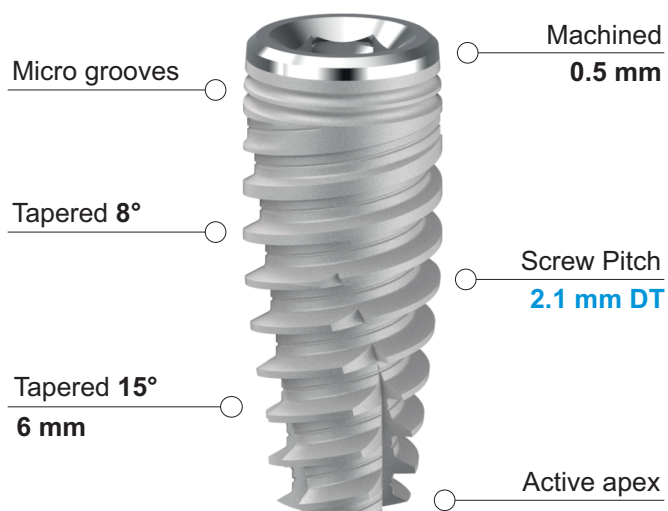
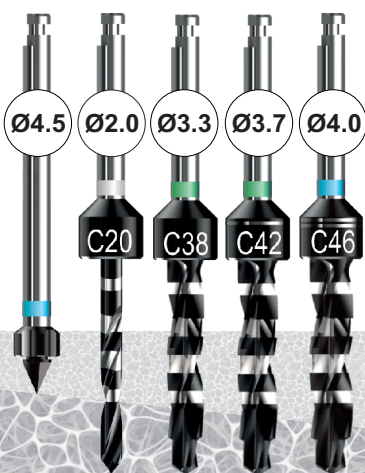
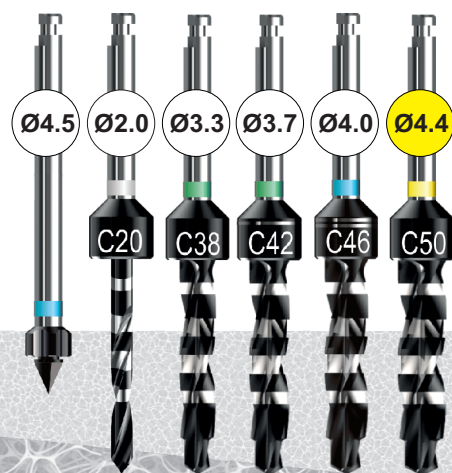
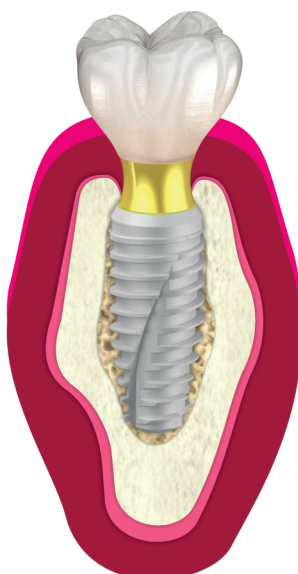
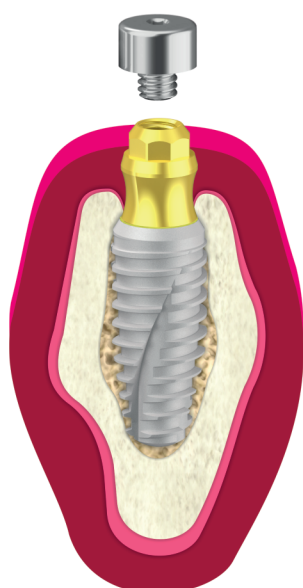
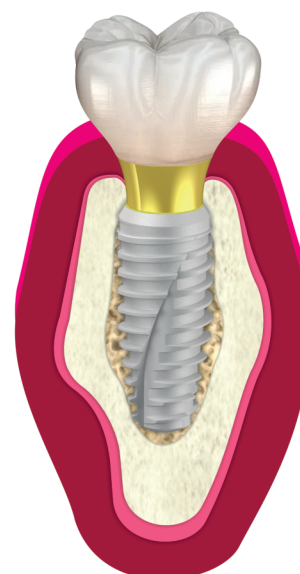
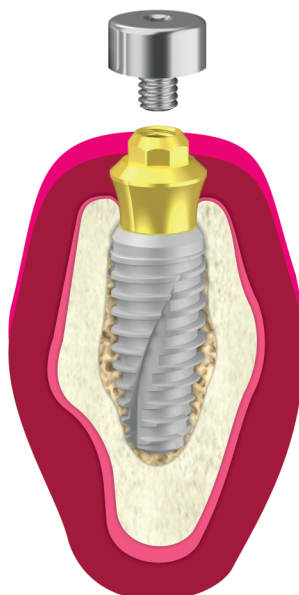
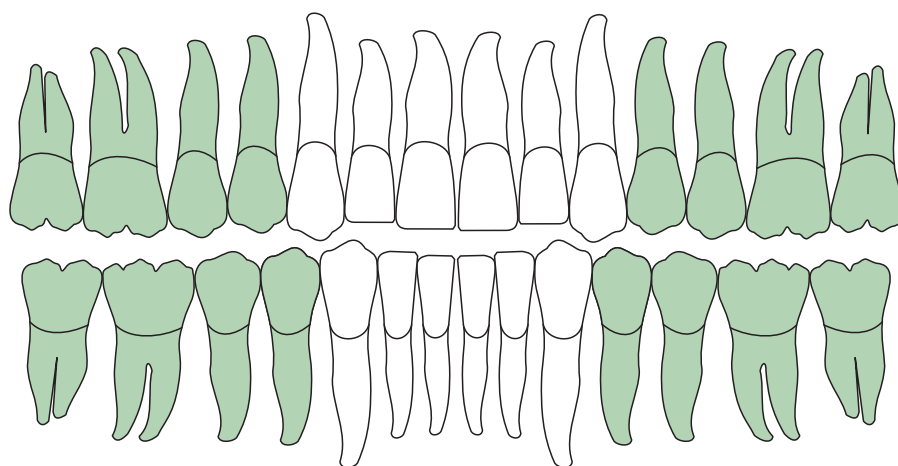
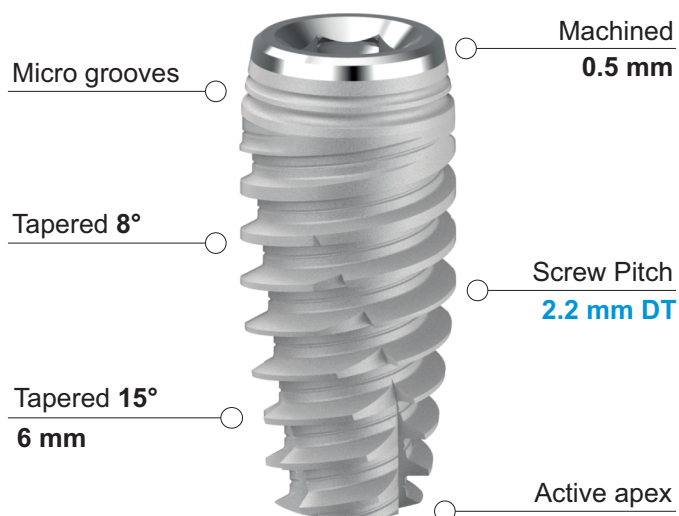
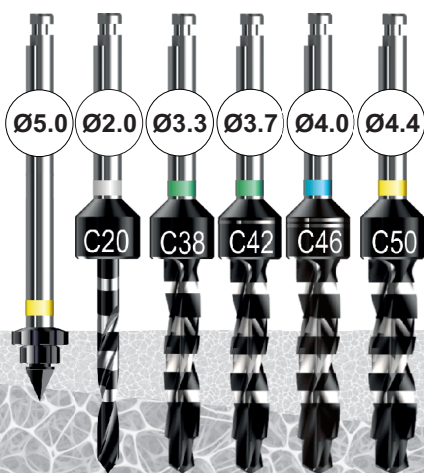
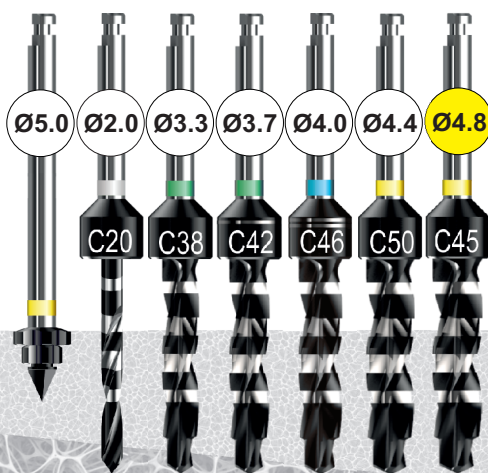
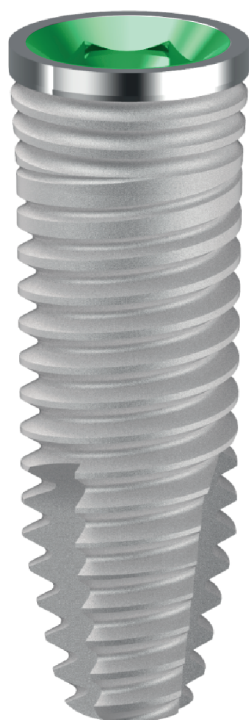
IC

IC 38


CYLINDRICAL SHAPE Internal Hexagon

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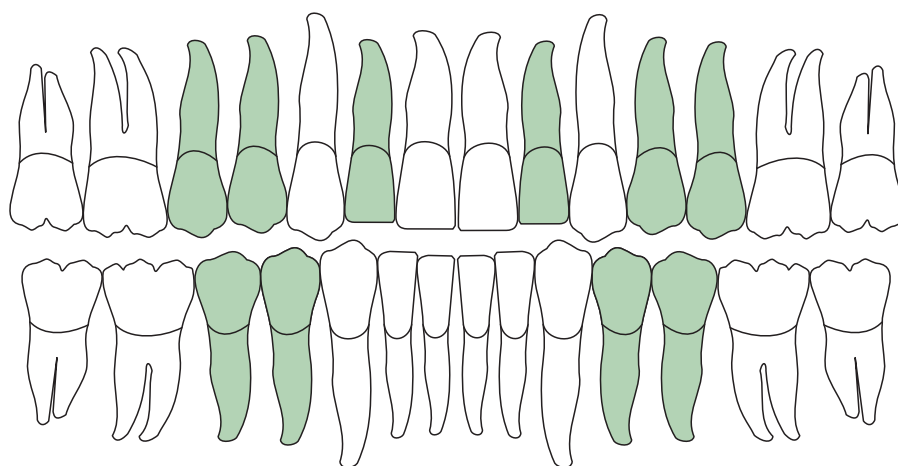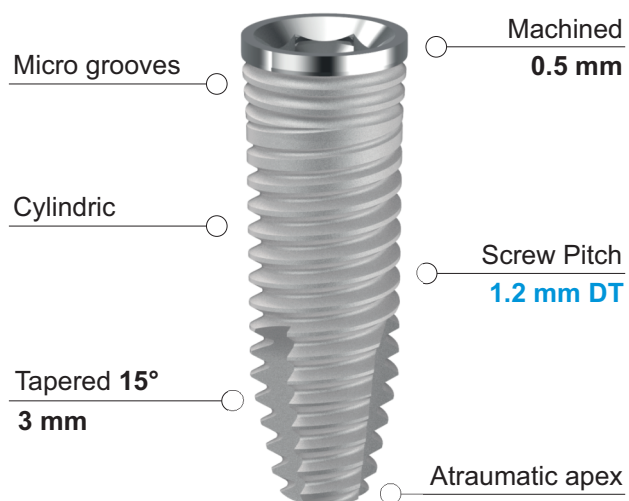
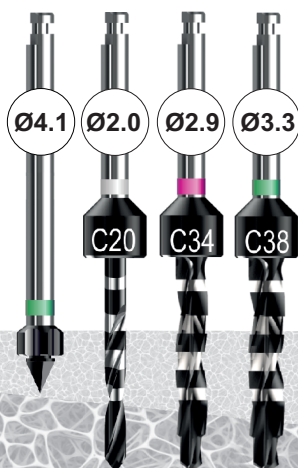
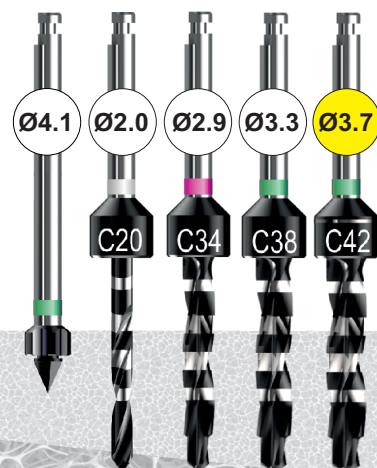
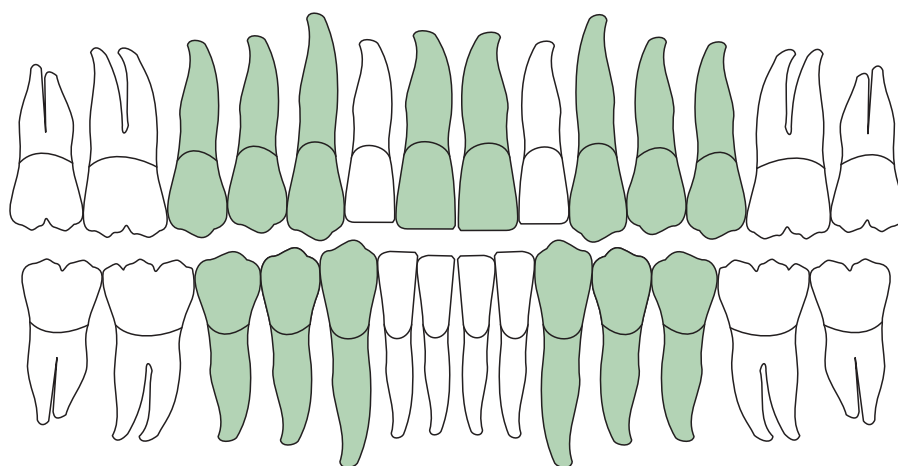
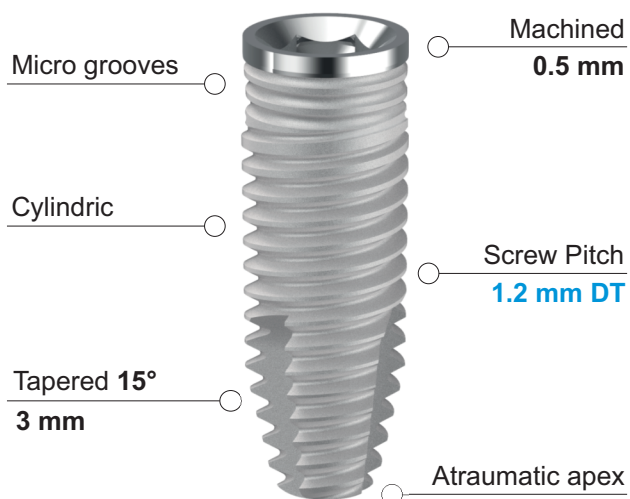
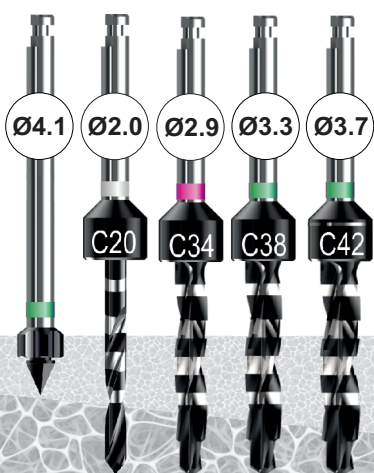
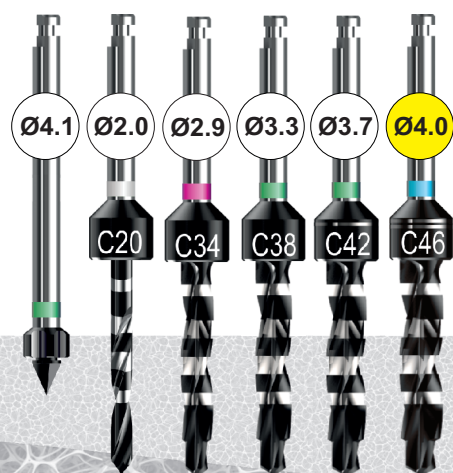
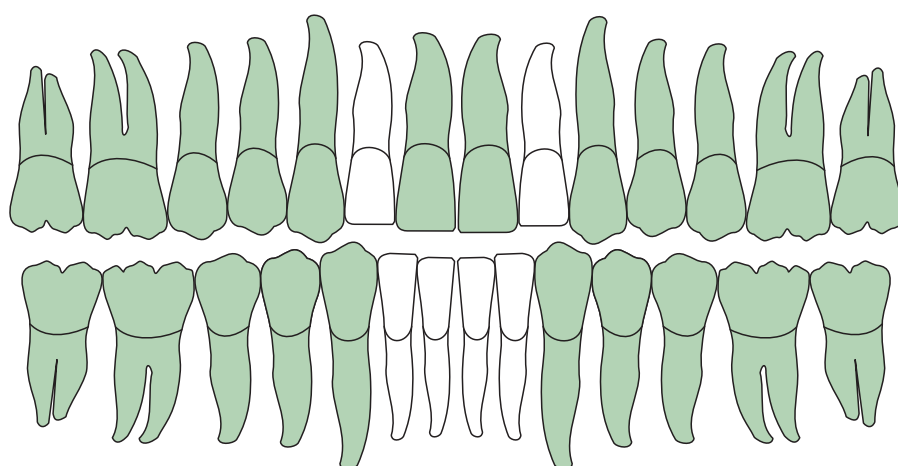
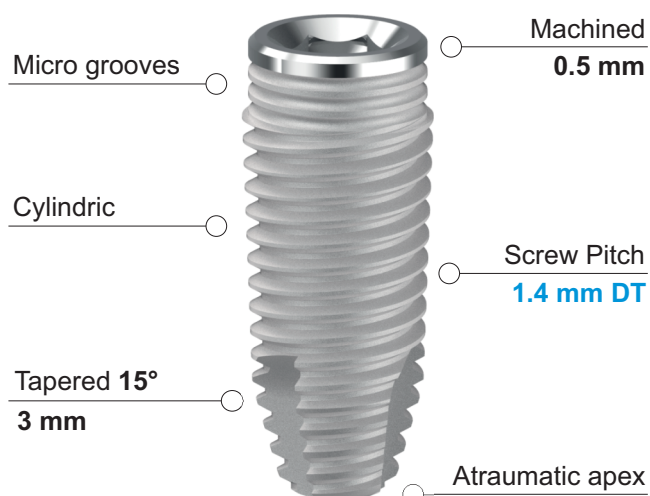
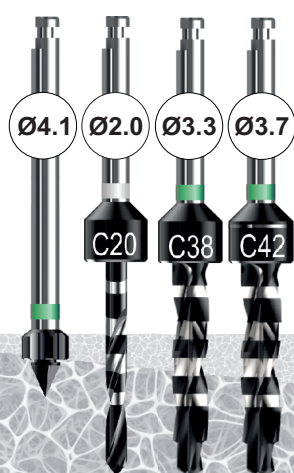
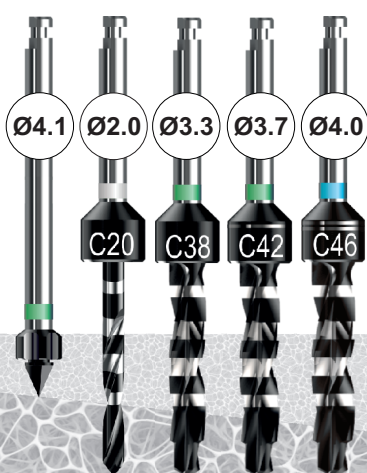
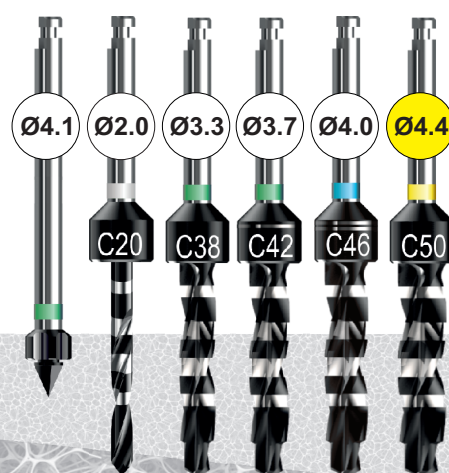
IC 50


LARGE CYLINDRICAL SHAPE Internal Hexagon

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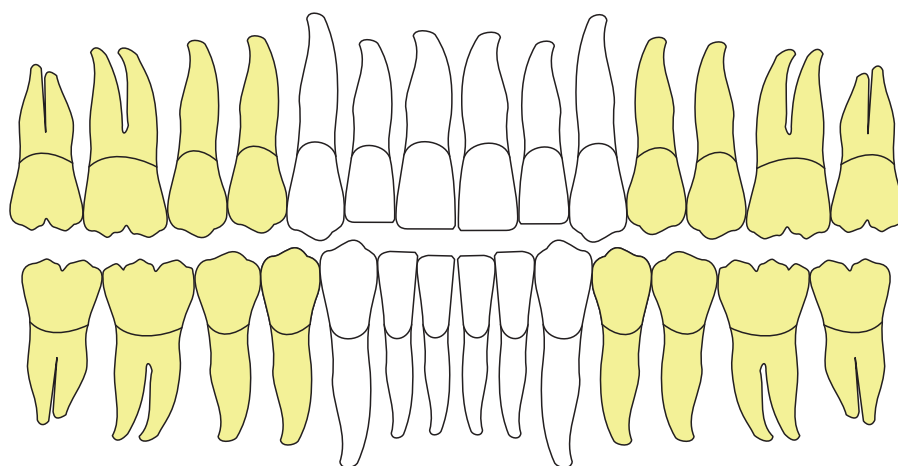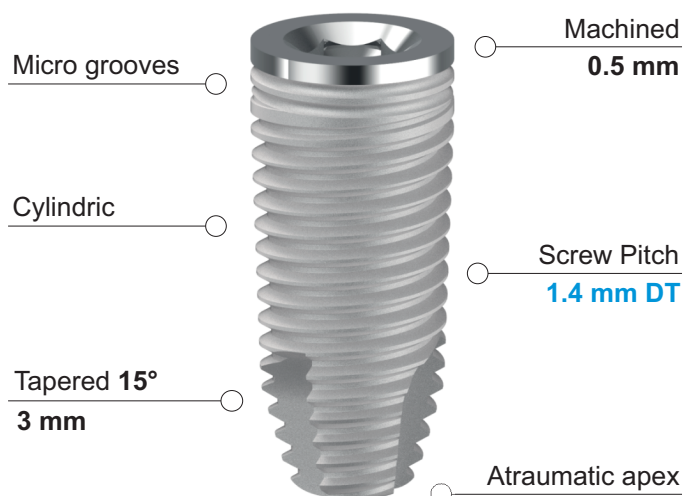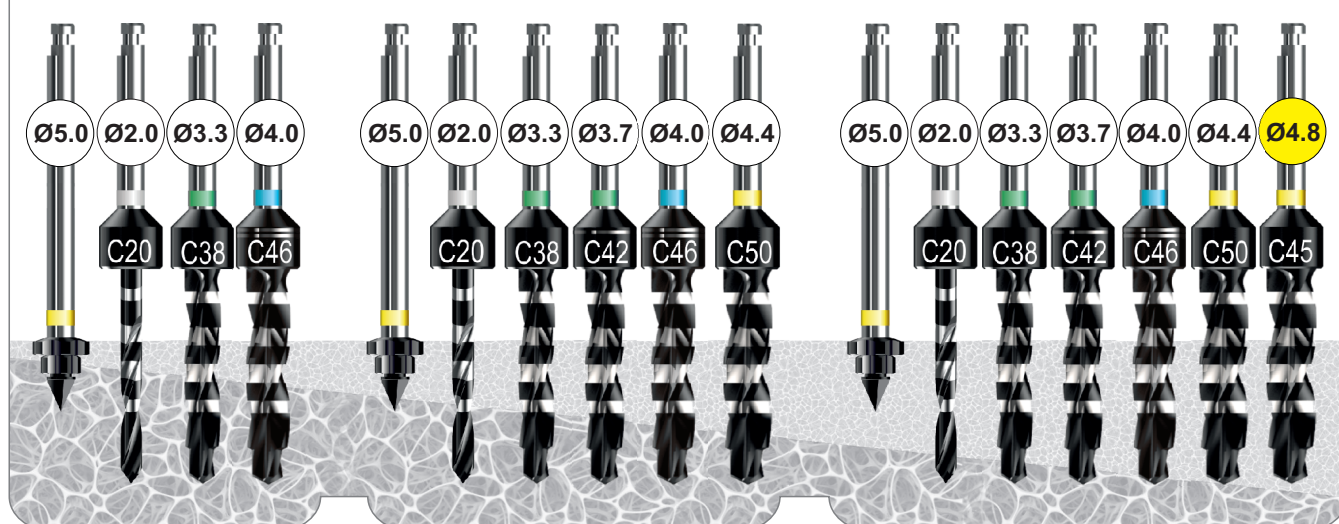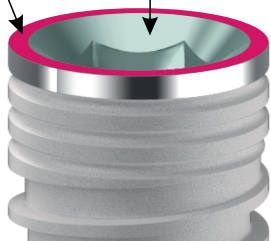
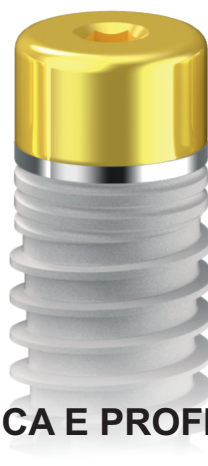
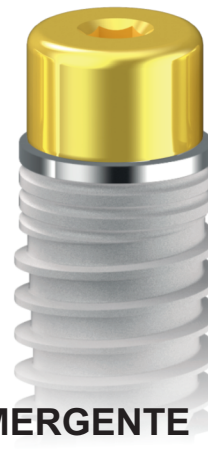
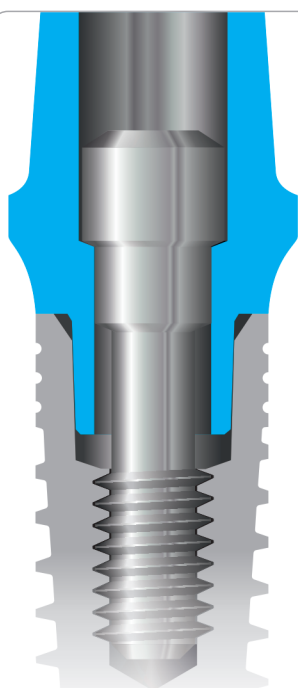
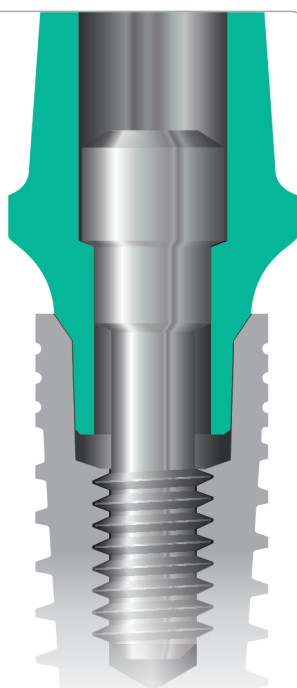
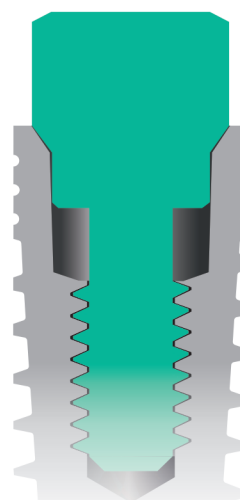
CPS Healing Abutment - Titanium Gr. 5 - Anodized

- IP 4HA2-W H 2 mm - **CPS** Standard Gingival Profile
- IP 4HA4-W H 4 mm - **CPS** Standard Gingival Profile
- IP 4HA6-W H 6 mm - **CPS** Standard Gingival Profile
- IP 4HA7-W H 7.5 mm - **CPS** Standard Gingival Profile

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FTF Healing Abutment - Titanium Gr. 5 - Anodized

- IP 4HA2 H 2 mm - **FTF** Standard Gingival Profile
- IP 4HA4 H 4 mm - **FTF** Standard Gingival Profile
- IP 4HA6 H 6 mm - **FTF** Standard Gingival Profile

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FTF Large Healing Abutment - Titanium Gr. 5 - Anodized

- IP 4HA2-E H 2 mm - **FTF** Large Gingival Profile
- IP 4HA4-E H 4 mm - **FTF** Large Gingival Profile
- IP 4HA6-E H 6 mm - **FTF** Large Gingival Profile

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


CPS Cover Screw - Titanium Gr. 4

- IP 4CS1 CPS Cover Screw - 1.5mm - Narrow Gingival Profile
- IP 4CS3 CPS Cover Screw - 3.0mm - Narrow Gingival Profile
- IP 4CS5 CPS Cover Screw - 5.0mm - Narrow Gingival Profile

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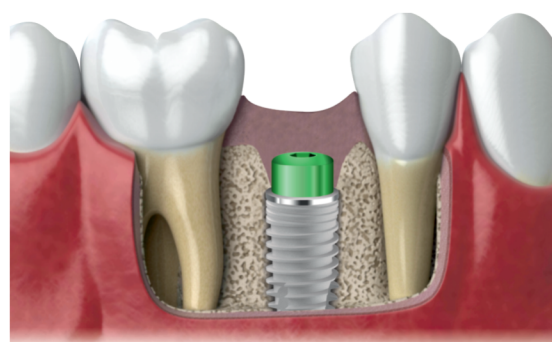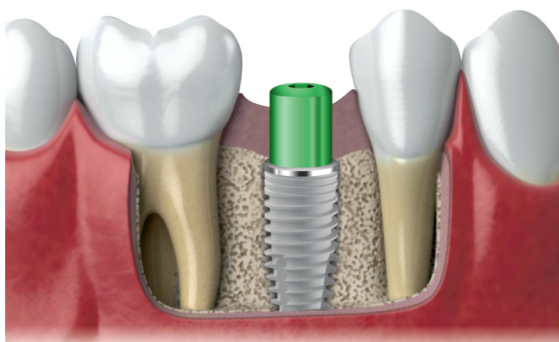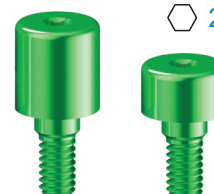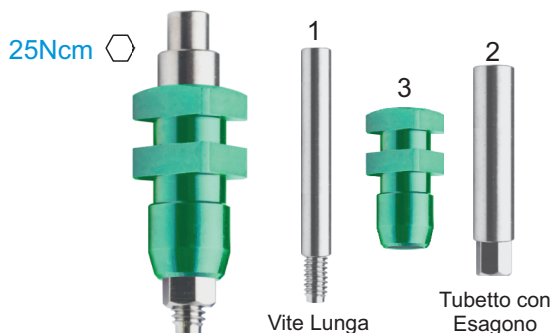
Transfer per impronta Pick Up - Titanio Gr. 5 **Pick Up Impression Coping - Titanium Gr. 5**

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Temporary Glueing Abutments - Titanium Gr. 5

IP 4TA-W Temporary Glueing Abutment - Antirotation Hex

IP 4TR-W Temporary Glueing Abutment - Rotating

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DM 1/2 Standard Hexagonal Driver Ø 1.27

Pilastro Temporaneo da Saldatura - Titanio Gr. 5

Temporary Welding Abutments - Titanium Gr. 5

IP 4TRS Temporary Welding Abutment - Rotating

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DM 1/2 Standard Hexagonal Driver Ø 1.27

Moncone Titanio Angolato - Titanio Gr. 5 - PVD TiN


Angled Titanium Abutment - Titanium Gr. 5 - PVD TiN

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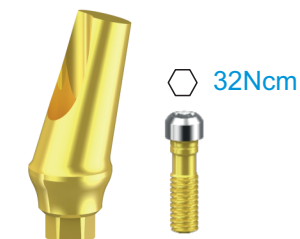
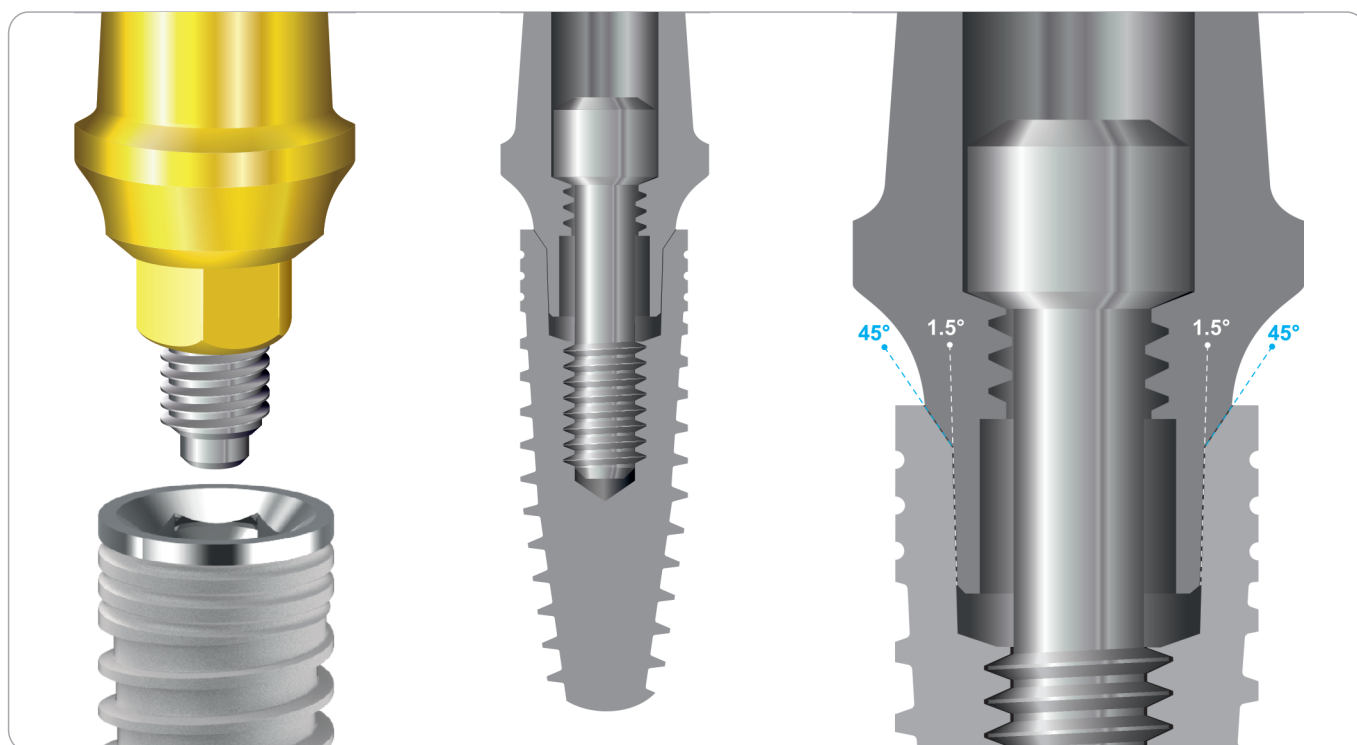
32Ncm

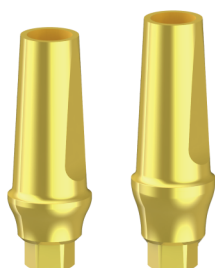Moncone Titanio Dritto - Titanio Gr. 5 - PVD TiN

Straight Titanium Abutment - Titanium Gr. 5 - PVD TiN

IP 4ST2-WF 0° - Straight H2mm - **CPS** Friction Fit
IP 4ST3-WF 0° - Straight H3.5mm - **CPS** Friction Fit
IP 5ST2-F 0° - Straight H2mm (Large) - **FTF** Friction Fit

IP VITG
DM 1/2

Vite di Ricambio Gold - Gold Spare Screw
 Standard Hexagonal Driver Ø 1.27

32Ncm

Moncone in Titanio Angolato 15° - Titanio Gr. 5 - PVD TiN

15° Angled Titanium Abutment - Titanium Gr. 5 - PVD TiN

IP 4AT152-WF 15° Angled H2mm - **CPS** Friction Fit
IP 4AT153-WF 15° Angled H3mm - **CPS** Friction Fit
IP 5ST152-F 15° - Angled H2mm (Large) - **FTF** Friction Fit

IP VITG
DM 1/2

Vite di Ricambio Gold - Gold Spare Screw
 Standard Hexagonal Driver Ø 1.27

32Ncm

Moncone in Titanio Angolato 25° - Titanio Gr. 5 - PVD TiN

25° Angled Titanium Abutment - Titanium Gr. 5 - PVD TiN

IP 4AT252-WF 25° Angled H2mm - **CPS** Friction Fit
IP 4AT253-WF 25° Angled H3mm - **CPS** Friction Fit
IP 5ST252-F 25° - Angled H2mm (Large) - **FTF** Friction Fit

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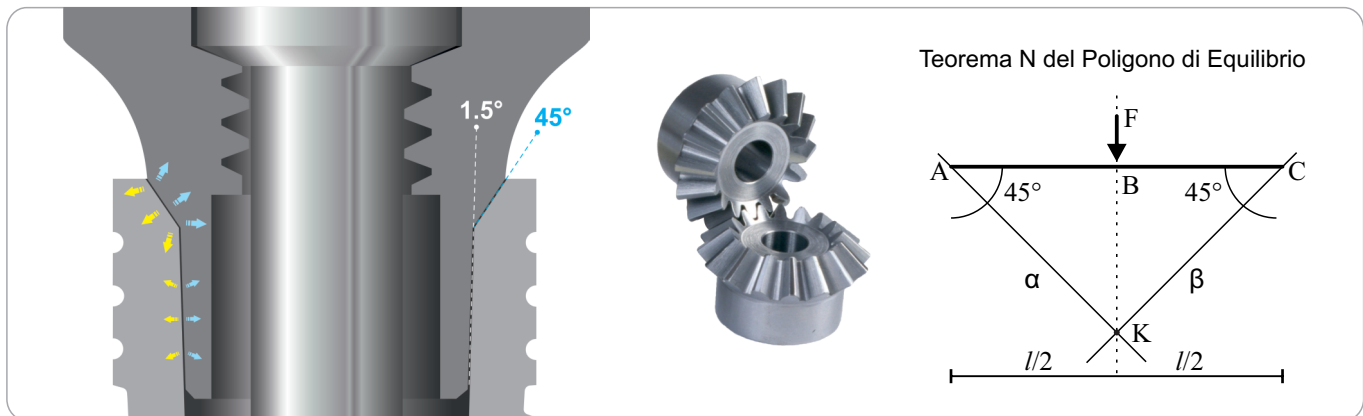
45° CONICAL SHOULDER CLOSURE

The 45° internal prosthetic support, already used in internal Zimmer type connections since the 90s, has an indisputable bibliography and moreover guarantees the distribution of vector forces that minimize the lateral bending moment.


This type of conical support, normally used in the distribution mechanics, increases the contact surface in relation to a plane surface and releases the forces on a curve of multi-directional vectors.

FRICTION FIT 1.5° CONICAL HEX CONNECTION

The 1.5° conicity on the flat surfaces of the male / female hexagon causes a so close reduction of the mechanical interface that confers a direct binding (known as "cold welding") between the abutment and the implant.

This effect is fully activated at the end of the prosthetic tightening of the passing screw, that reaches 30 Ncm of torque and guarantees a bacterial seal through the reduction of the interface gap to zero.

In order to remove the abutment it will be necessary to use an alternative removal screw, that needs to be screwed in the original one place: its unique function is to push the abutment from the inside vertically and to remove it from its bound position.

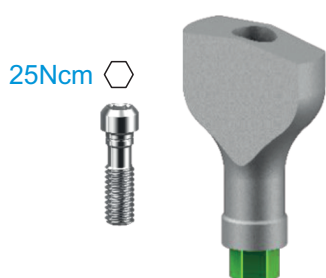
Transfer per Impronta Digitale - Titanio Gr. 5 + MS
Scan Abutment - Titanium Gr. 5 + Micro Sandblasting
IP 4SA2 DS Antirotation Hex - **Narrow** Gingival Profile

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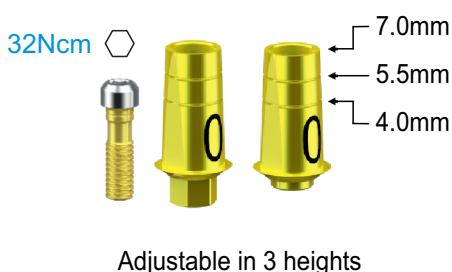Flat Digital Link for Bonding Technique - Titanium Gr. 5

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Open Links for Gluing for Angled Screwing - Titanium Gr. 5

IP4TBAI0DS-W Gingival Profile **H 0.6mm** - Antirotation Hex

IP4TBAI1DS-W Gingival Profile **H 1.5mm** - Antirotation Hex

IP4TBAI3DS-W Gingival Profile **H 3.0mm** - Antirotation Hex


IP4TBRAI0DS-W Gingival Profile **H 0.6mm** - Rotating

IP4TBRAI1DS-W Gingival Profile **H 1.5mm** - Rotating

IP4TBRAI3DS-W Gingival Profile **H 3.0mm** - Rotating

MDT6DS Manual Torque Driver for angled hole

DCT6DS RA Contr-angle Torque Driver for angled hole

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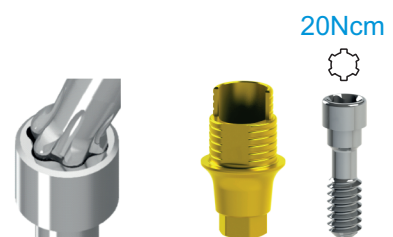Analogical Link for Bonding Technique - Titanium Gr. 5


IP 4SZB Analogical Titanium Link - Antirotation Hex - **CPS**

IP VIT Vite di Ricambio - Spare Screw

DM 1/2 Standard Hexagonal Driver Ø 1.27

oneabutment octa prosthetics

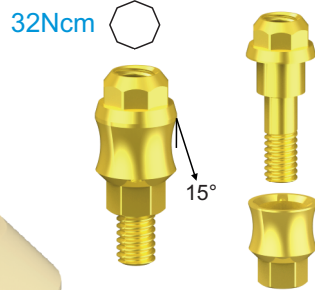
Oneabutment Octa 3.8 - Titanio Gr. 5

Oneabutment Octa 3.8 - Titanium Gr. 5

IP4O 381	Oneabutment Transgengivale H1.5 mm
IP4O 382	Oneabutment Transgengivale H2.5 mm
IP4O 383	Oneabutment Transgengivale H3.5 mm
IP4O 384	Oneabutment Transgengivale H4.5 mm
IP4O 385	Oneabutment Transgengivale H5.5 mm

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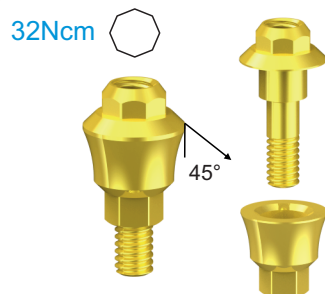Oneabutment Octa 4.8 - Titanio Gr. 5

Oneabutment Octa 4.8 - Titanium Gr. 5

IP4O 481	Oneabutment Transgengivale H1.5 mm
IP4O 482	Oneabutment Transgengivale H2.5 mm
IP4O 483	Oneabutment Transgengivale H3.5 mm
IP4O 484	Oneabutment Transgengivale H4.5 mm
IP4O 485	Oneabutment Transgengivale H5.5 mm

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Standard Hexagonal Driver Ø 1.20

Transfer universale per impronta Pick Up - Titanio Gr. 5

Universal Pick Up Impression Coping - Titanium Gr. 5

EPIAO 38	Impression Coping Octa Ø 3.8 mm e Ø 4.8 mm
----------	--

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EPBALO 48DS Analogo Digitale Octa Ø 4.8 mm

Pilastro Temporaneo Octa - Titanio Gr. 5


Octa Temporary Abutment - Titanium Gr. 5

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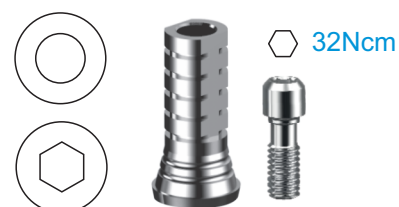
Short Bonding Technique Octa Base Link - Titanium Gr. 5

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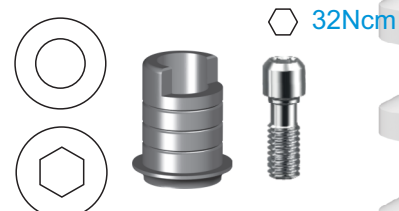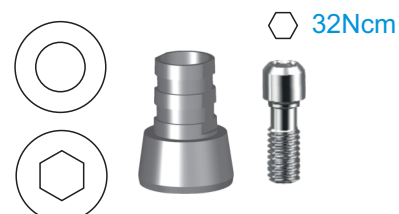
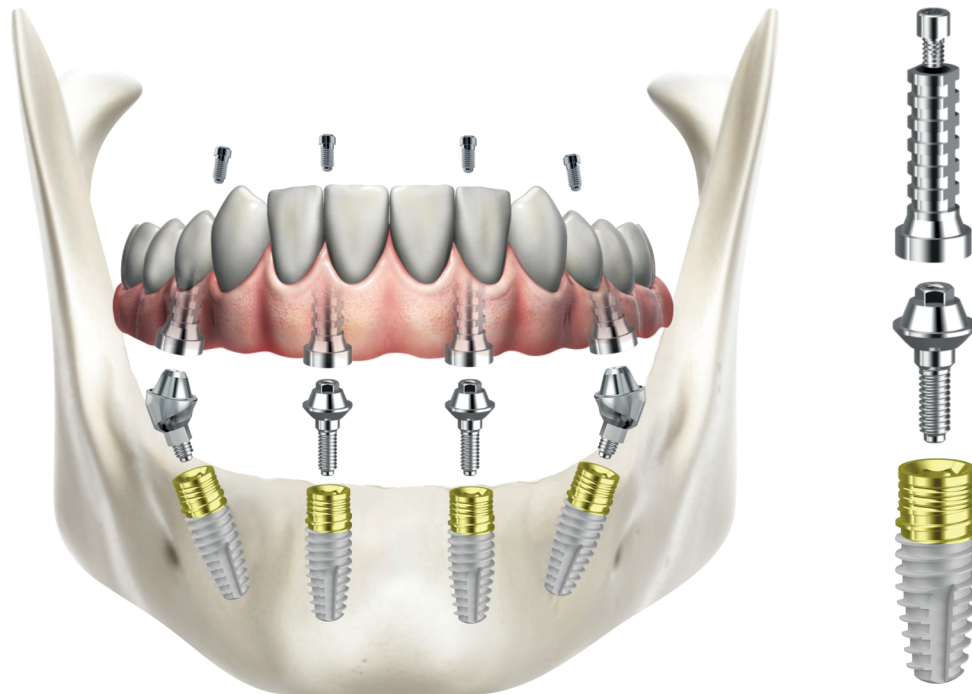
32Ncm

Pilastro Toronto Angolato 30° - Titanio Gr. 5

30° Toronto Angled Abutment - Titanium Gr. 5

- IP MA303-W 30° - H 3 mm - **CPS** Standard Gingival Profile
- IP MA304-W 30° - H 4 mm - **CPS** Standard Gingival Profile
- IP MA305-W 30° - H 5 mm - **CPS** Standard Gingival Profile

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toronto universal prosthetics

Pilastro Temporaneo Toronto - Titanio Gr. 5

Toronto Temporary Abutment - Titanium Gr. 5

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DMX 1/2 Standard Torx Driver

Pilastro Calcinabile Toronto - Acrilico

Toronto Castable Abutment - Acrylic

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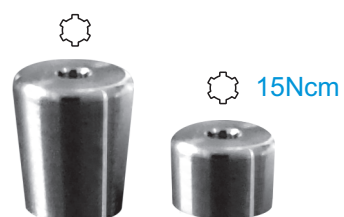
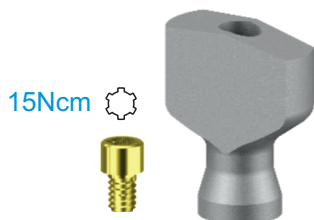
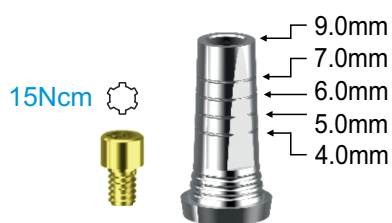
EP MAL2DS H 10 mm - Universal Shape - Whit screw

EP MAL2DS-V H 10 mm - Universal Shape - Whitout screw

DM 1/2 Standard Hexagonal Driver Ø 1.27

Viti di Ricambio Toronto - Titanio Gr. 5 **Toronto Spare Screws - Titanium Gr. 5**

EP MV Prosthetic Spare Screw
EP MVL Impression Spare Screw

DMX 1/2 Standard Torx Driver

universal prosthetics

Pilastro a Sfera per Overdenture - Titanio Gr. 5 - PVD TiN

Overdenture Ball Abutment - Titanium Gr. 5 - PVD TiN

IP 4BA1-W	H 1 mm - CPS Standard Gingival Profile - NORMO
IP 4BA2-W	H 2 mm - CPS Standard Gingival Profile - NORMO
IP 4BA3-W	H 3 mm - CPS Standard Gingival Profile - NORMO
IP 4BA4-W	H 4 mm - CPS Standard Gingival Profile - NORMO
IP 4BA6-W	H 6 mm - CPS Standard Gingival Profile - NORMO


DM 1/2 Standard Hexagonal Driver Ø 1.27

OT Equator per Overdenture - Titanio Gr. 5 - PVD TiN

Overdenture OT Equator - Titanium Gr. 5 - PVD TiN

IP 4EQM1	H 1 mm - Only Abutment - Narrow Gingival Profile
IP 4EQM2	H 2 mm - Only Abutment - Narrow Gingival Profile
IP 4EQM3	H 3 mm - Only Abutment - Narrow Gingival Profile
IP 4EQM4	H 4 mm - Only Abutment - Narrow Gingival Profile
IP 4EQM5	H 5 mm - Only Abutment - Narrow Gingival Profile
IP 4EQM6	H 6 mm - Only Abutment - Narrow Gingival Profile

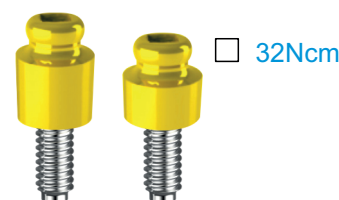DM EQ Squared Equator Driver

Kit OT Equator per Overdenture - Titanio Gr. 5 - PVD TiN

Overdenture OT Equator Kit - Titanium Gr. 5 - PVD TiN

IP 4EQ1	H 1 mm - Narrow Gingival Profile
IP 4EQ2	H 2 mm - Narrow Gingival Profile
IP 4EQ3	H 3 mm - Narrow Gingival Profile
IP 4EQ4	H 4 mm - Narrow Gingival Profile
IP 4EQ5	H 5 mm - Narrow Gingival Profile
IP 4EQ6	H 6 mm - Narrow Gingival Profile


DM EQ Squared Equator Driver

EasyFlag3D advanced digital impression

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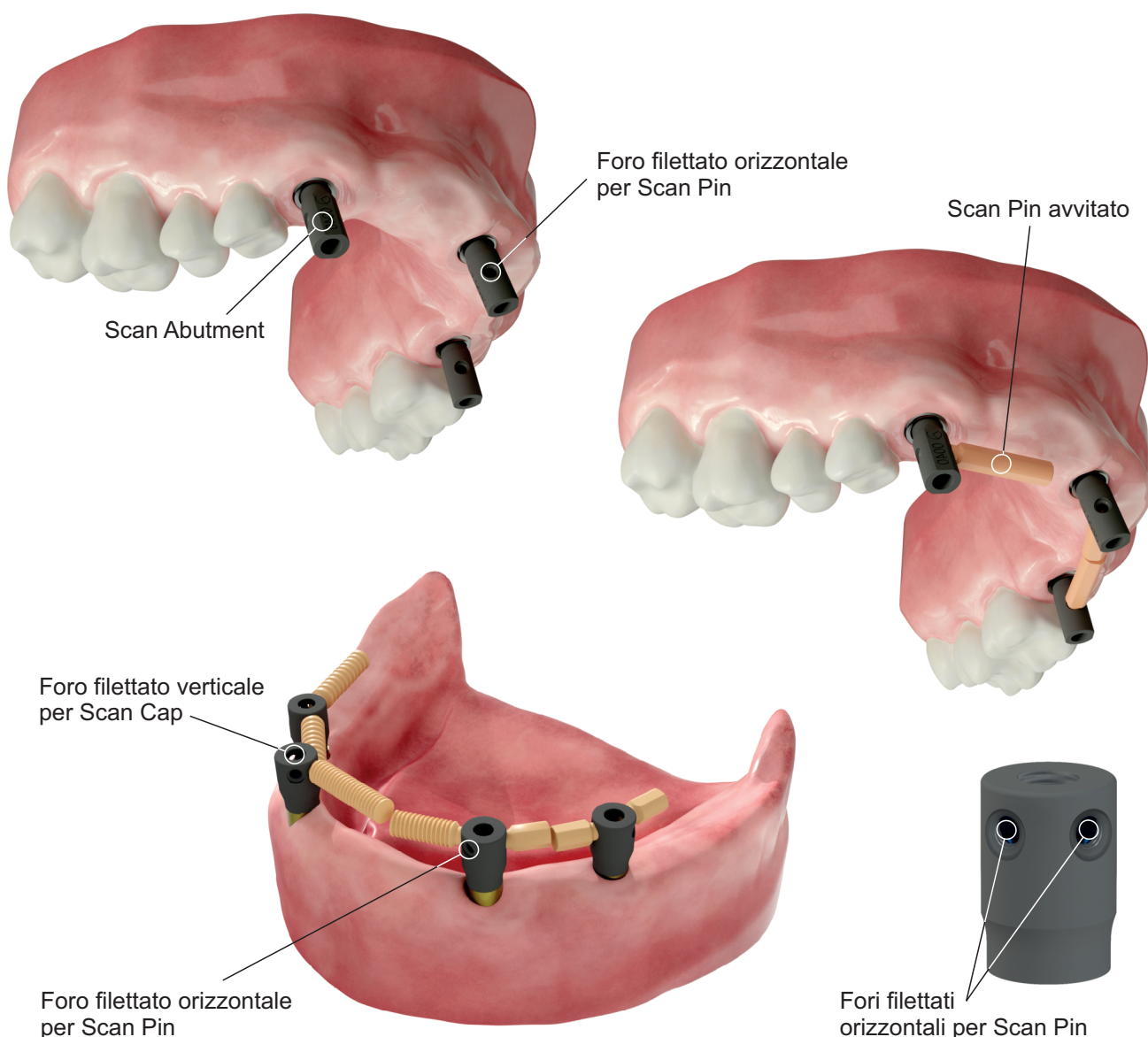The main limitation in digital impression precision was due to the distance between the scanabutments.

The lack of a continuous profile, combined with the smooth gingival tissue, required the software to continually re-evaluate the position, resulting in frequent errors and system crashes that compromised the final accuracy.


EasyFlag3D is a customizable continuous intraoral scanning system that ensures fast STL image acquisition and accuracy.

EasyFlag3D advanced digital impression

EF3D Scan Abutment - Titanio Gr. 5

Scan Abutment EF3D - Titanium Gr. 5


IP4SADS-EF	Scan Abutment EF3D Internal Hexagon CPS
EPMSADS-EF	Scan Abutment EF3D Toronto MUA
EPMSAQDS-EF	Scan Abutment EF3D UNI-Q-MUA

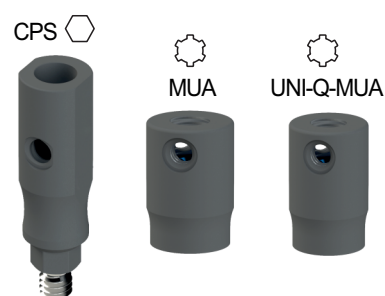EF3D Scan Pin Type A - Triangolare - PEEK Medicale

Type A Scan Pin EF3D - Triangular Shape - Medical PEEK


SPA6-EF	Scan Pin EF3D Type A 6 mm (5 pz)
SPA9-EF	Scan Pin EF3D Type A 9 mm (5 pz)
SPA13-EF	Scan Pin EF3D Type A 13 mm (5 pz)

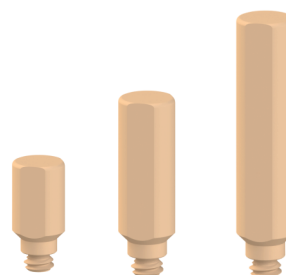EF3D Scan Pin Type B - Triangolare - PEEK Medicale

Type B Scan Pin EF3D - Triangular Shape - Medical PEEK


SPB6-EF	Scan Pin EF3D Type B 6 mm (5 pz)
SPB9-EF	Scan Pin EF3D Type B 9 mm (5 pz)
SPB13-EF	Scan Pin EF3D Type B 13 mm (5 pz)

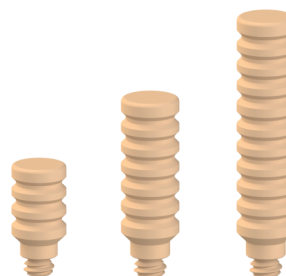EF3D Scan Pin Type C - Triangolare - PEEK Medicale

Type C Scan Pin EF3D - Triangular Shape - Medical PEEK

SPC6-EF	Scan Pin EF3D Type C 6 mm (5 pz)
SPC9-EF	Scan Pin EF3D Type C 9 mm (5 pz)
SPC13-EF	Scan Pin EF3D Type C 13 mm (5 pz)

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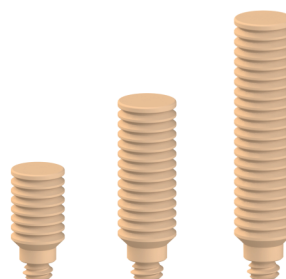Scan Caps EF3D - Occlusal Registration Screwcaps - Medical PEEK

SCR38-EF	Scan Cap EF3D Regular 3.8 mm
SCR60-EF	Scan Cap EF3D Regular 6 mm
SCR80-EF	Scan Cap EF3D Regular 8 mm

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Healing Abutment - Titanium Gr. 5 - Anodized

IP 5HA3 H 3 mm - Large Gingival Profile

IP 5HA5 H 5 mm - Large Gingival Profile

DM 1/2

Standard Hexagonal Driver Ø 1.27

32Ncm 

Moncone Calcinabile - Acrilico

Castable Abutment - Acrylic

IP 5CA Antirotation Hex - Large Gingival Profile

IP 5CR Rotating - Large Gingival Profile

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32Ncm 

Moncone in Titanio Dritto - Titanio Gr. 5 - PVD TiN

Straight Titanium Abutments - Titanium Gr. 5 - PVD TiN

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32Ncm 

Monconi in Titanio Angolati - Titanio Gr. 5 - PVD TiN

Angled Titanium Abutments - Titanium Gr. 5 - PVD TiN

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Ucla Gold Abutment - Gold/Platinum/Iridium Base

IP 5GA Antirotation Hex - Large Gingival Profile

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Pilastrì Temporanei - Titanio Gr. 5

Temporary Abutments - Titanium Gr. 5

IP 4TA-W Glueing Temporary Abutment - Antirotation Hex

IP 4TR-W Glueing Temporary Abutment - Rotating

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DM 1/2 Standard Hexagonal Driver Ø 1.27

Pilastrì Temporanei da Saldatura - Titanio Gr. 5

Welding Temporary Abutments - Titanium Gr. 5

IP 4TRS Welding Temporary Abutment - Rotating

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
Aesthetic Temporary Abutment - Medical PEEK

IP 5PA Antirotation Hex - Large Gingival Profile

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DM 1/2 Standard Hexagonal Driver Ø 1.27

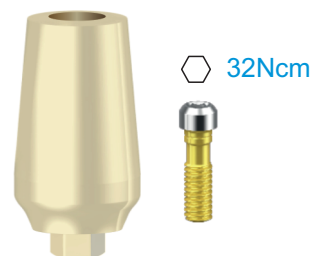Transfer per impronta Pick Up - Titanio Gr. 5


Pick Up Impression Coping - Titanium Gr. 5

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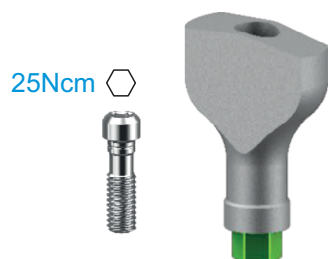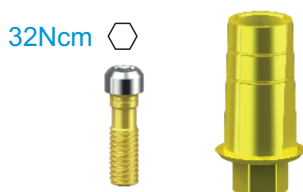
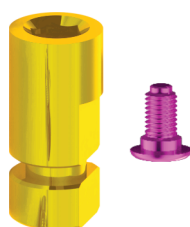
Laboratory Digital Analog - Titanium Gr. 5 - Anodized

IP 5AL2DS H 12 mm - Universal Shape

Pilastro a Sfera per Overdenture - Titanio Gr. 5 - PVD TiN

Overdenture Ball Abutment - Titanium Gr. 5 - PVD TiN

IP 4BA1-W H 1 mm - **CPS** Standard Gingival Profile - NORMO

IP 4BA2-W H 2 mm - **CPS** Standard Gingival Profile - NORMO

IP 4BA3-W H 3 mm - **CPS** Standard Gingival Profile - NORMO

IP 4BA4-W H 4 mm - **CPS** Standard Gingival Profile - NORMO

IP 4BA6-W H 6 mm - **CPS** Standard Gingival Profile - NORMO

DM 1/2 Standard Hexagonal Driver

Kit OT Equator per Overdenture - Titanio Gr. 5 - PVD TiN

Overdenture OT Equator Kit - Titanium Gr. 5 - PVD TiN

IP 4EQ1 H 1 mm - Narrow Gingival Profile

IP 4EQ2 H 2 mm - Narrow Gingival Profile

IP 4EQ3 H 3 mm - Narrow Gingival Profile

IP 4EQ4 H 4 mm - Narrow Gingival Profile

IP 4EQ5 H 5 mm - Narrow Gingival Profile

DMEQ Squared Equator Driver

toronto universal prosthetics

Pilastro Toronto Dritto - Titanio Gr. 5

Toronto Straight Abutment - Titanium Gr. 5

IP MS1-W H 1 mm - Standard Gingival Profile

IP MS2-W H 2 mm - Standard Gingival Profile

IP MS3-W H 3 mm - Standard Gingival Profile

DCDM Toronto Straight Abutment Driver

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Toronto Healing Abutment and Impression Coping

EP MI Rotating - Standard Gingival Profile

EP MHA H 4 mm - Standard Gingival Profile

EP MHA6 H 6 mm - Large Gingival Profile

DMX 1/2 Standard Torx Driver

SPARE SCREWS & MOUNTERS



REPLACEMENT PRODUCTS

SCREWS & MOUNTERS


spare screws

Esagono Interno IK / IA / IC - Internal Hexagon IK / IA / IC

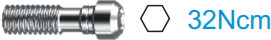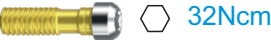	IP VIT20	H 20mm Transfer Screw - Ø 1.8mm
	IP VIT30	H 30mm Transfer Screw - Ø 1.8mm
	IP VIT	H 8mm Prosthetic Screw - Ø 1.8mm
	IP VITG	Prosthetic Screw - Ø 1.8mm - PVD TiN
	IP VITT6DS	Prosthetic Screw T6 Avvitamento Ang. Ø 1.8mm
	IP VITMA	Angled 17°/30° MUA Abutment Screw Ø 1.8mm
	EP MV	Toronto Prosthetic Spare Screw - Ø 1.4mm
	EP MVL	Toronto Impression Spare Screw - Ø 1.4mm
	EPVO	Vite Protesica OCTA Prosthetic Spare Screw

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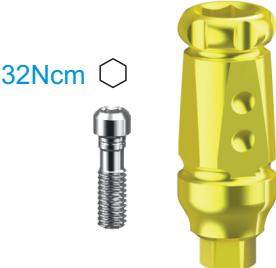		MTA Moulder Transfer Abutment - Titanio Gr. 5	
		MTA Moulder Transfer Abutment - Titanium Gr. 5	
		IPM 40	H 12mm - Ø 4.7mm Internal Hex. (34 - 38 - 42)
		IPM 46	H 12mm - Ø 5.5mm Internal Hex. (46 - 50)
		DM 1/2	Standard Hexagonal Driver

INSTRUMENTS & DRILLS

IMPLANT TOOLS

⑤ standard surgical drills

FRESE HBM - Hard Bone Management

	07001604	Cortical Drill - Ø 1.6 mm - 11 mm
	04002104	Cortical Drill - Ø 2.1 mm - 11 mm
	FLS2313	Lance Drill - Ø 2.3 mm - 13 mm

FIS- Corticals

	FIS 30	Cortical Drill - Ø 3.0 mm
	FIS 35	Cortical Drill - Ø 3.7 mm
	FIS 40	Cortical Drill - Ø 4.1 mm
	FIS 45	Cortical Drill - Ø 4.5 mm
	FIS 50	Cortical Drill - Ø 5.0 mm

FCS - Cylindricals

	FCS 20	Cylindrical Drill - Ø 2.0 mm - Long 18 mm
	FCS 30	Cylindrical Drill - Ø 2.6 mm - Long 18 mm
	FCS 34	Cylindrical Drill - Ø 2.9 mm - Long 18 mm
	FCS 38	Cylindrical Drill - Ø 3.3 mm - Long 18 mm
	FCS 42	Cylindrical Drill - Ø 3.7 mm - Long 18 mm
	FCS 46	Cylindrical Drill - Ø 4.0 mm - Long 18 mm
	FCS 50	Cylindrical Drill - Ø 4.4 mm - Long 18 mm
	FCS 55	Cylindrical Drill - Ø 4.8 mm - Long 18 mm

⑤ standard surgical drills

FKS - Tapered

FKS 30	Tapered Drill - Ø 2.4 mm - Long 18 mm
FKS 34	Tapered Drill - Ø 2.9 mm - Long 18 mm
FKS 38	Tapered Drill - Ø 3.3 mm - Long 18 mm
FKS 42	Tapered Drill - Ø 3.7 mm - Long 18 mm
FKS 46	Tapered Drill - Ø 4.0 mm - Long 18 mm
FKS 50	Tapered Drill - Ø 4.4 mm - Long 18 mm
FKS 55	Tapered Drill - Ø 4.8 mm - Long 18 mm
FKS 60	Tapered Drill - Ø 5.2 mm - 15 mm
FKS 80	Tapered Drill - Ø 7.2 mm - 15 mm


F- Nasal Pterygoid Long Drills

FCC 2025	N. Pterygoids Standard Drill - Ø 2.0mm
FKC 3425	N. Pterygoids Standard Drill - Ø 2.9mm
FKC 3825	N. Pterygoids Standard Drill - Ø 3.3mm
FKC 4225	N. Pterygoids Standard Drill - Ø 3.6mm

Mucotomi - Tissue Punches

ST M13	Manual Tissue Punche - Ø3.0mm
ST M14	Manual Tissue Punche - Ø4.1mm
ST M15	Manual Tissue Punche - Ø5.1mm
ST M23	Contrangle Tissue Punche - Ø3.0mm
ST M24	Contrangle Tissue Punche - Ø4.1mm
ST M25	Contrangle Tissue Punche - Ø5.1mm

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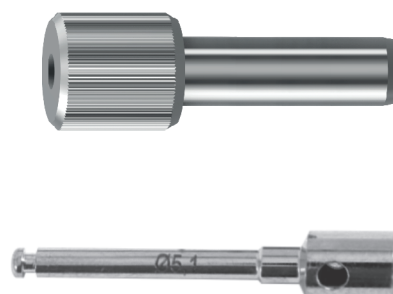FBP-OT6 Bone Profiler OneTime - Ø 6.0 mm

PF Drill Extension - Prolunga per Frese

FG - Guided Nasal Pterygoid Long Drills

FGR 2025L N. Pterygoids Guided Drill - Ø 2.0mm

FGR 3425L N. Pterygoids Guided Drill - Ø 2.9mm

FGR 3825L N. Pterygoids Guided Drill - Ø 3.3mm

FGR 4225L N. Pterygoids Guided Drill - Ø 3.6mm

COMPUTER GUIDED SURGERY

GUIDED SURGERY TOOLS


computer guided surgery

MASCHERINA CHIRURGICA - Mascherina on demand multi-impianti

SURGICAL GUIDE - Multi-implants on demand surgical guide

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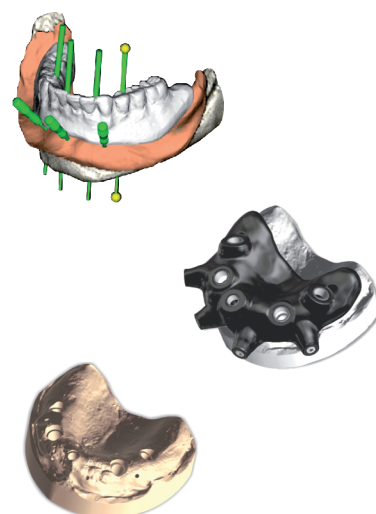
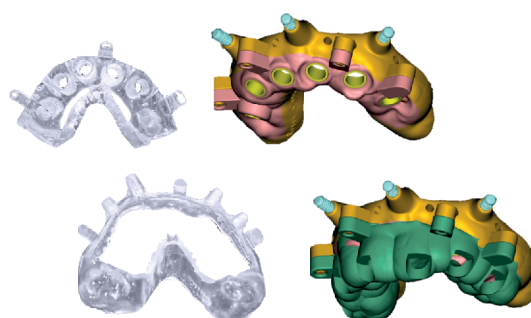
MODELLI - Modelli Stereolitografati

MODELS - Stereolithographed models

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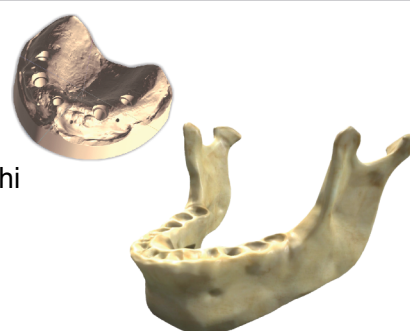
IMMEDIATE TEMPORARY PROSTHESIS - Drilled PMMA and Milled Titanium Framework

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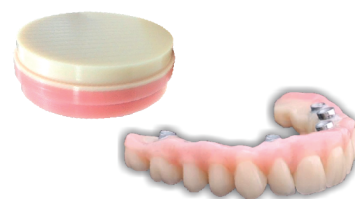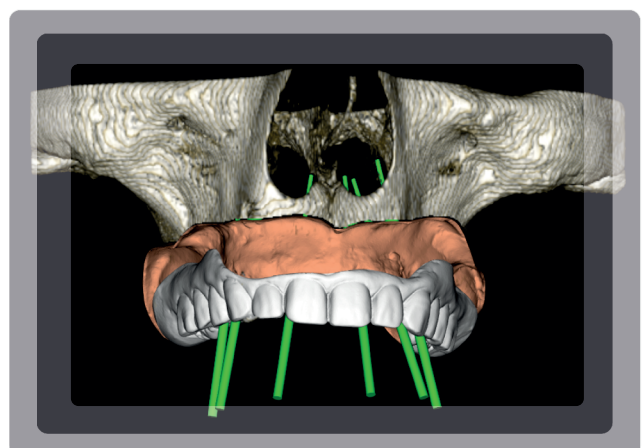
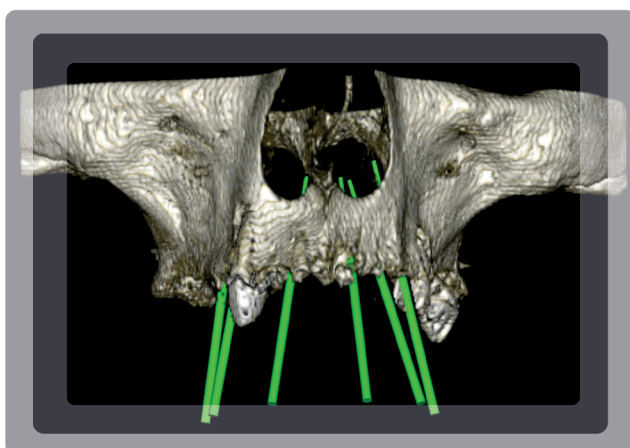
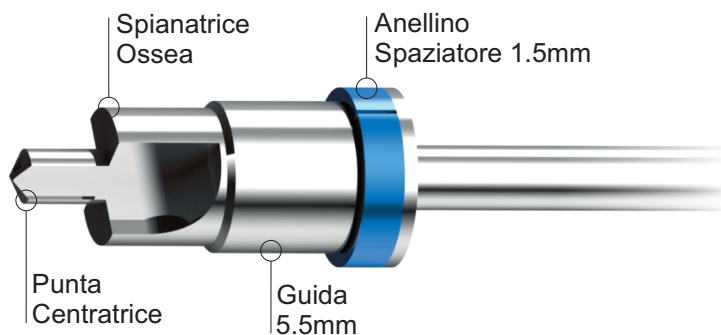
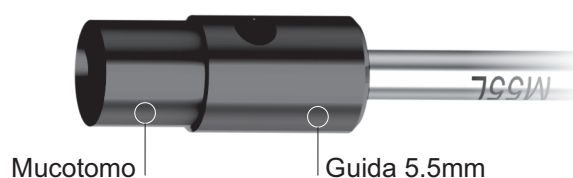
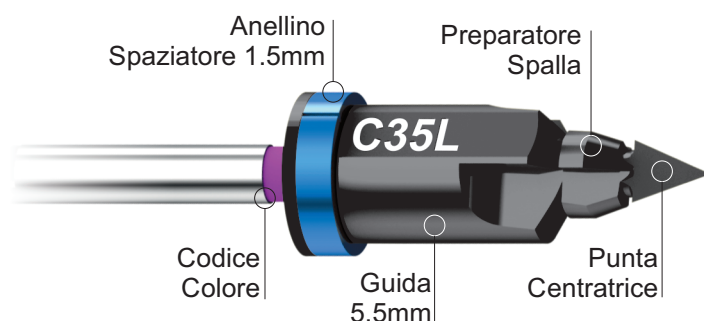
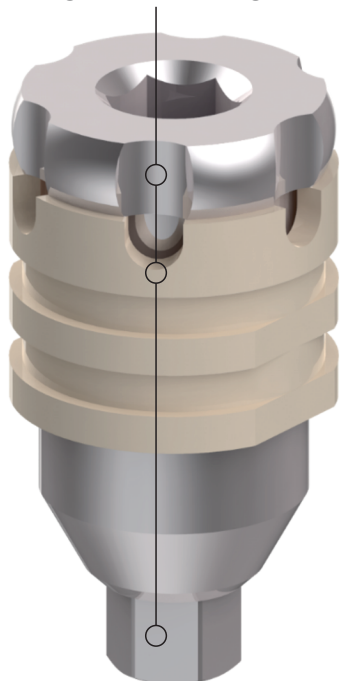
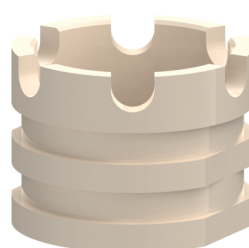
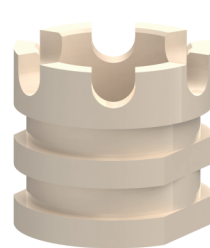
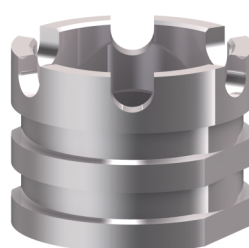
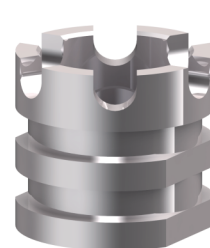
Guided Mounters and Sleeves

IP4MG 4S

Mounter Guidata - Internal Hex ● Small Drive

IP4MG 4L

Mounter Guidata - Internal Hex ● Large Drive

BOC41P

Boccola 4.15 Small **PEEK** - Sleeve 4.15 Small **PEEK**

BOC41T

Boccola 4.15 Small **Ti** - Sleeve 4.15 Small **Ti**

BOC55P

Boccola 5.5 Large **PEEK** - Sleeve 5.5 Large **PEEK**

BOC55T

Boccola 5.5 Large **Ti** - Sleeve 5.5 Large **Ti**

BOC15

Boccola 4.15 Small **PEEK** - Sleeve 4.15 Small **PEEK**

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
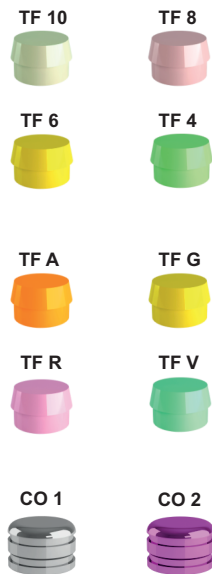
Kit LOCATOR® Internal Hexagon - Abutment + Box + Caps

IP4LC1	Kit LOCATOR h1 mm - Ø4.0 CPS	●
IP4LC2	Kit LOCATOR h2 mm - Ø4.0 CPS	●
IP4LC3	Kit LOCATOR h3 mm - Ø4.0 CPS	●
IP4LC4	Kit LOCATOR h4 mm - Ø4.0 CPS	●
IP4LC5	Kit LOCATOR h5 mm - Ø4.0 CPS	●
IP4LC6	Kit LOCATOR h6 mm - Ø4.0 CPS	●

LOCATOR®

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LOCATOR® Accessories and Spare Parts

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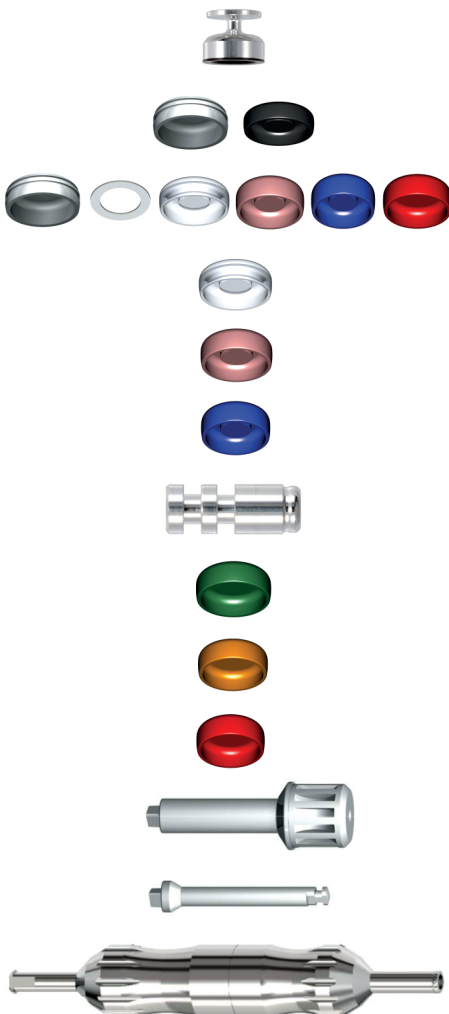
instruments & tools

Cacciaviti Manuali - Manual Screw Drivers

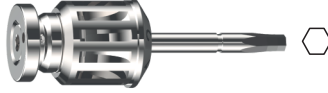	DM 01	Hexagonal 1.27 - Short - 12mm
	DM 02	Hexagonal 1.27 - Medium - 17mm
	DM 03	Hexagonal 1.27 - Long - 27mm
	DMX 01	Torx - Short - 15mm
	DMX 02	Torx - Long - 20mm
	DMRP120S	Hexagonal 1.20 - Short
	DMRP120M	Hexagonal 1.20 - Medium
	DMRP120L	Hexagonal 1.20 - Long
	DM T6DS	T6 Angled Torx Driver - 15mm
	DMEQ	Equator Squared - Short - 15mm

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	DC 01	Hexagonal 1.27 - Short - 15mm
	DC 02	Hexagonal 1.27 - Medium - 20mm
	DC 03	Hexagonal 1.27 - Long - 30mm
	ACDRCA	Hexagonal 1.20
	DC T6DS	T6 Angled Torx Driver - 20mm
	DC T6DSXL	T6 Angled Torx Driver - 40mm
	DCEQ	Equator Squared - Medium - 20mm
	DCX 01	Torx - Short - 15mm
	DCX 02	Torx - Long - 20mm


Cacciaviti Pilastri Transgengivali - OCTA Transgingival Abutment Drivers

	DCM	Straight Toronto MUA Contra-angle Driver
	MU10AM	Straight Toronto MUA Manual Driver
	DCOC01	Octagonal One Abutment Driver - Short
	DCOC02	Octagonal One Abutment Driver - Long
	ACOC01	Octagonal One Abutment Driver - Short
	ACOC02	Octagonal One Abutment Driver - Long

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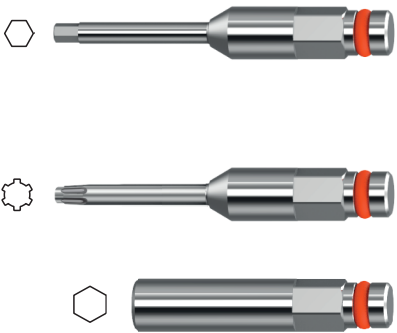
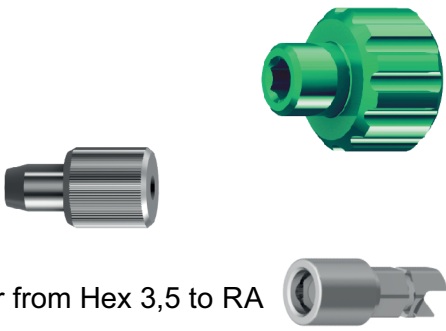
Prolunghe da Cricchetto IK/IC/IA/EC - IK/IC/IA/EC Ratchet Implant Drivers

- PC 01R** Short - Hexagonal Implant Driver
- PC 02R** Medium - Hexagonal Implant Driver
- PC 03R** Long - Hexagonal Implant Driver
- PC 02** Medium - Hexagonal Implant Driver **No O-Ring**
- PC 03** Long - Hexagonal Implant Driver **No O-Ring**
- PR PC** Prolunga per PC - PC Extension

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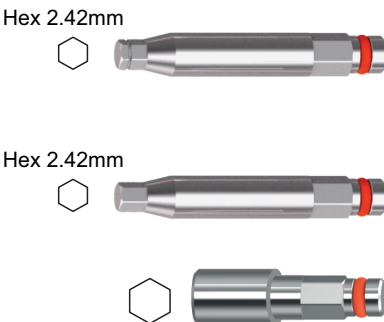
- AC 01** Short - Hexagonal Implant Driver
- AC 02** Medium - Hexagonal Implant Driver
- AC 03** Long - Hexagonal Implant Driver

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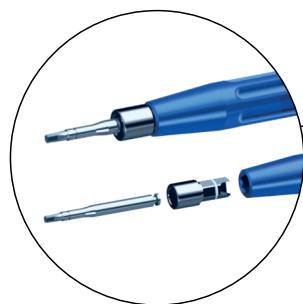Accessories for Manual Application

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
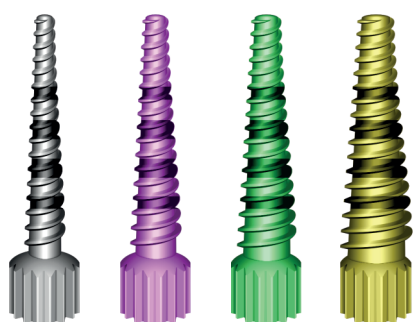
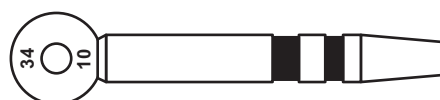
Bone Expanders and Compactors - Titanium Gr. 5

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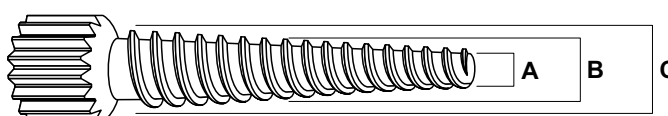
CRD2-G

CRD2 Cricchetto Dinamometrico 10/70 Ncm - Torque Controll

CRD2-G Ghiera Cricchetto CRD2 con Attacco PC - CRD2 Ratchet Wheel with PC Connection

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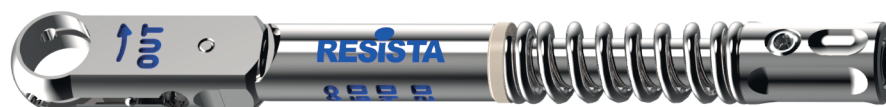
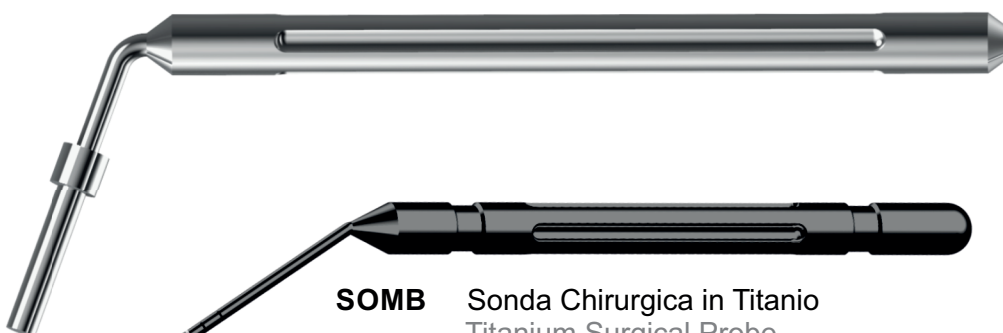
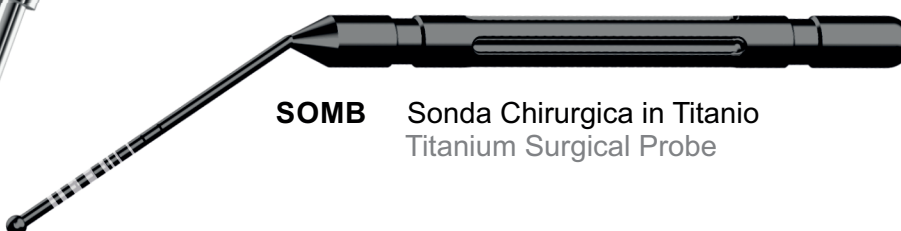
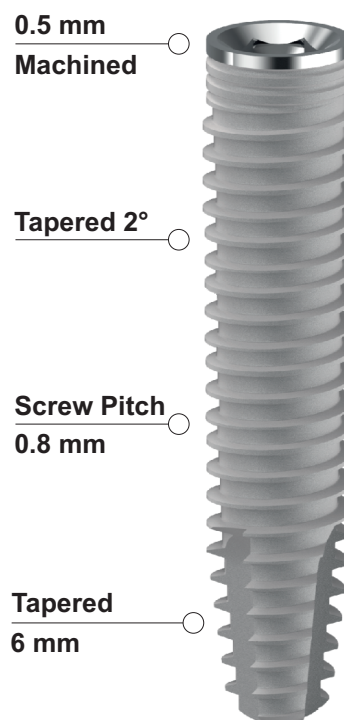
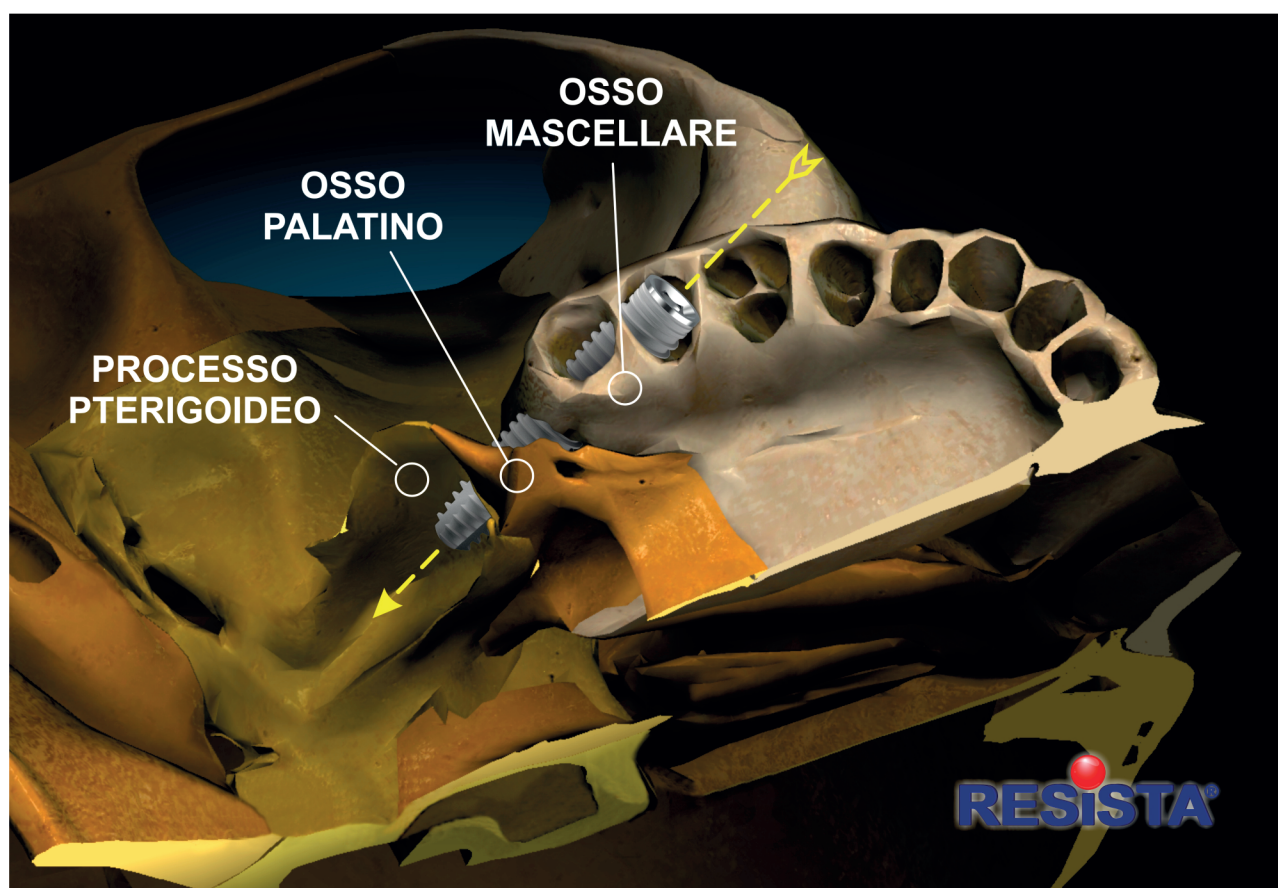
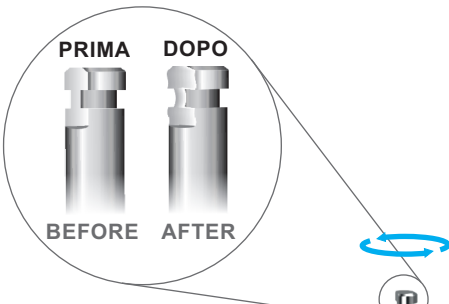
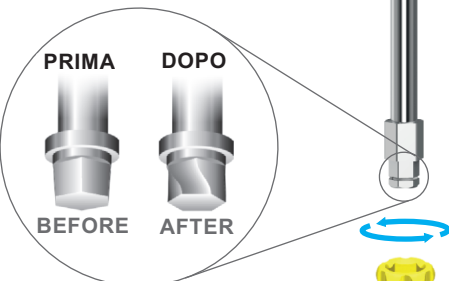
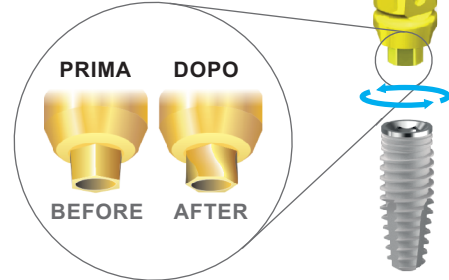
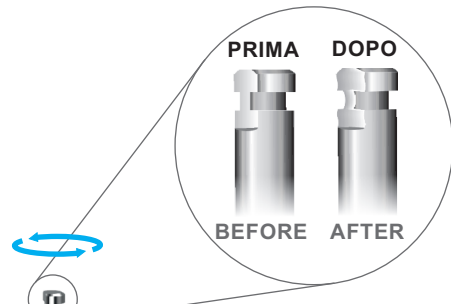
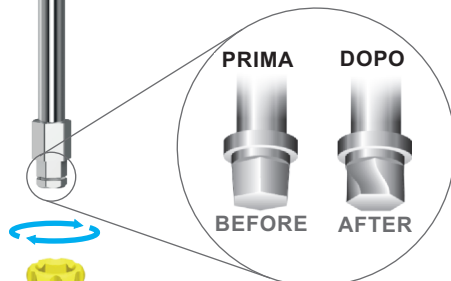
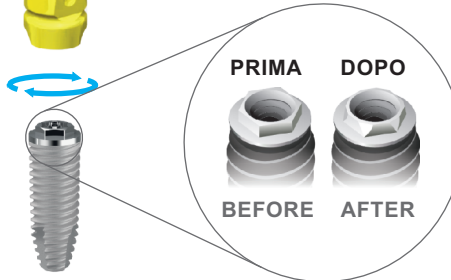
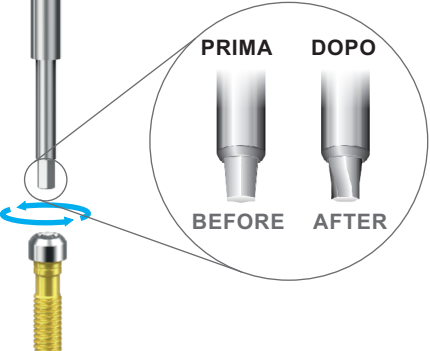
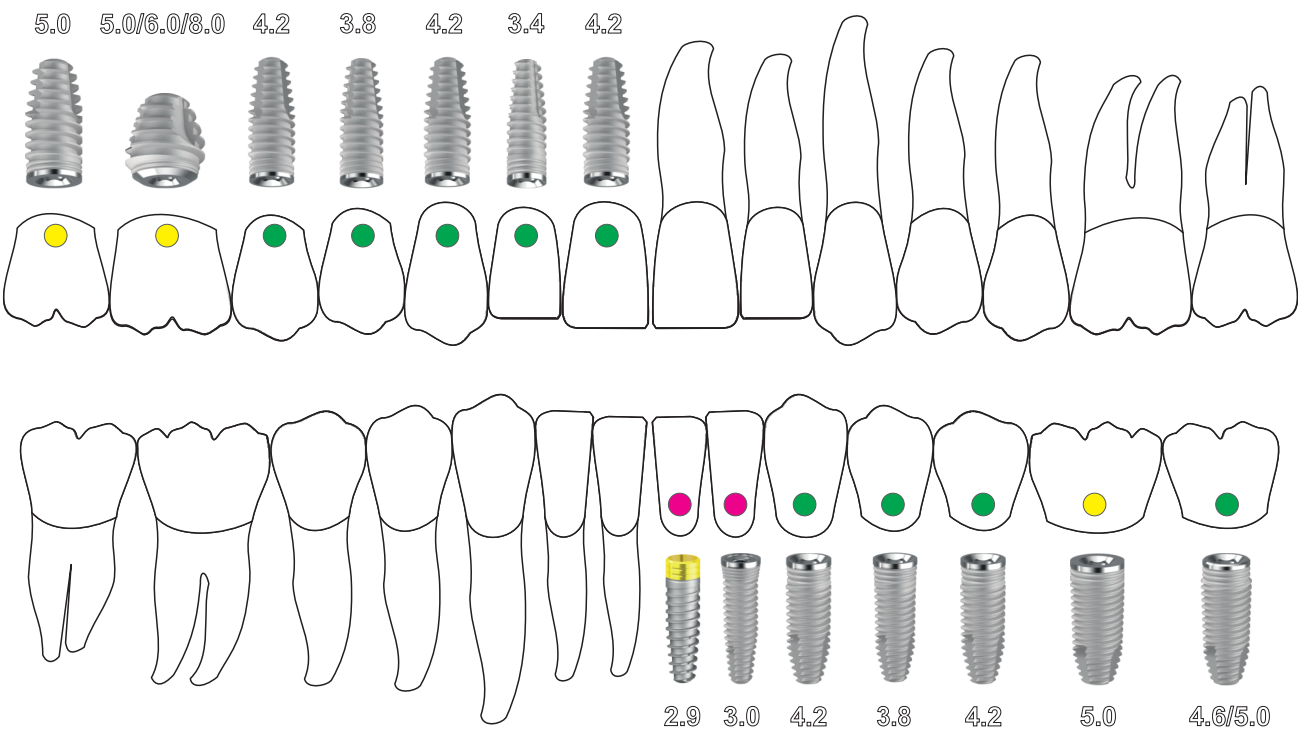
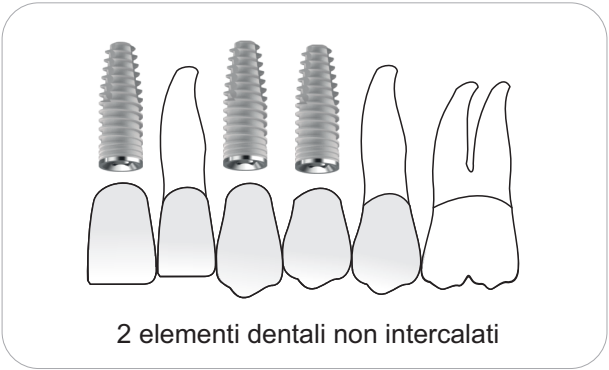
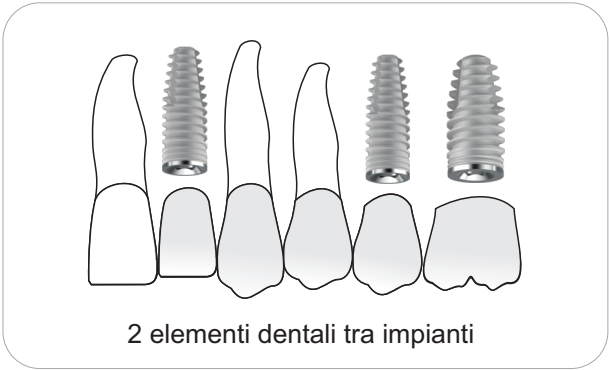
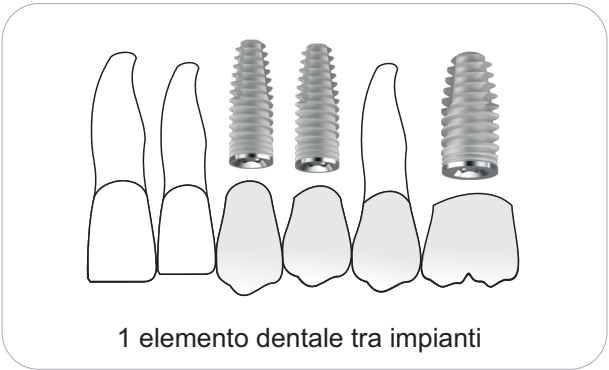
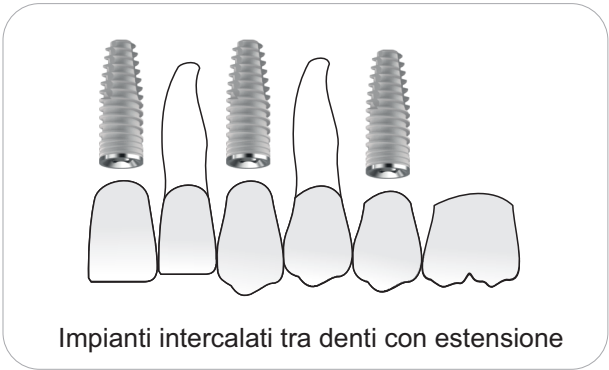
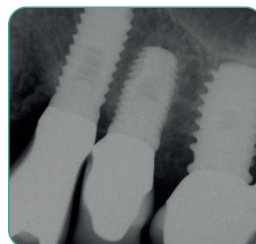
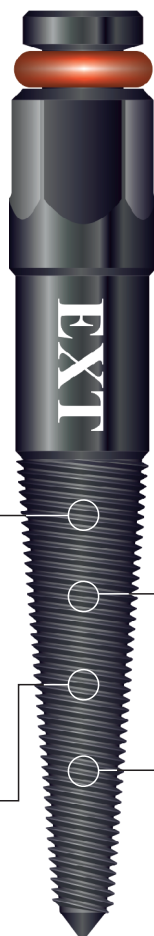
KIT 007 Implant Extractor Kit - 2 Extractors + 1 Box + CRND

Siringhe per Anestesia - Inox

Anesthesia Syringes - Stainless Steel

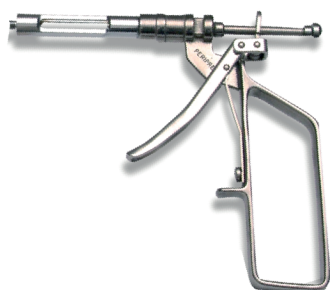PERIPRESS

22 PERIPRESS - Siringa Inox 1.8cc - Set 2 Portafiala

32 PERIPRESS - Siringa Standard 1.8cc - Set 2 Portafiala

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© calibrated surgical drills